



US011484660B2

(12) **United States Patent**  
**Creaturo**

(10) **Patent No.:** **US 11,484,660 B2**

(45) **Date of Patent:** **Nov. 1, 2022**

(54) **FORCE ACTUATED INJECTION DEVICE**

(71) Applicant: **Parenteral Technologies, LLC**, Siesta Key, FL (US)

(72) Inventor: **Michael A. Creaturo**, Siesta Key, FL (US)

(73) Assignee: **Parenteral Technologies, LLC**, Siesta Key, FL (US)

(\*) Notice: Subject to any disclaimer, the term of this patent is extended or adjusted under 35 U.S.C. 154(b) by 283 days.

(21) Appl. No.: **16/853,378**

(22) Filed: **Apr. 20, 2020**

(65) **Prior Publication Data**

US 2020/0246552 A1 Aug. 6, 2020

**Related U.S. Application Data**

(63) Continuation of application No. 15/200,298, filed on Jul. 1, 2016, now Pat. No. 10,625,026, which is a (Continued)

(51) **Int. Cl.**  
**A61M 5/315** (2006.01)  
**A61M 5/46** (2006.01)  
(Continued)

(52) **U.S. Cl.**  
CPC ..... **A61M 5/3158** (2013.01); **A61M 5/31511** (2013.01); **A61M 5/31555** (2013.01);  
(Continued)

(58) **Field of Classification Search**

CPC ..... A61M 5/46; A61M 2005/208; A61M 5/3243; A61M 2005/3267; A61M 5/3158;  
(Continued)

(56) **References Cited**

U.S. PATENT DOCUMENTS

4,927,416 A 5/1990 Tomkiel  
5,267,972 A 12/1993 Anderson  
(Continued)

OTHER PUBLICATIONS

International Search Report dated Apr. 9, 2015 for International Application No. PCT/US15/12152.

*Primary Examiner* — Bhisma Mehta


*Assistant Examiner* — John J Crawford

(74) *Attorney, Agent, or Firm* — Suiter Swantz pc llo

(57) **ABSTRACT**

A force actuated injection device including a body adapted to contain medication, a grip slidably disposed over the body, a lock disposed between the body and the grip, the lock when locked preventing movement of the grip relative to the body and when unlocked permitting movement of the grip relative to the body, a needle guard retractable relative to the body, a biasing member biasing the needle guard away from the body, and a resistance band disposed between the body and the needle guard adapted to resist initial needle guard retraction relative to the body, wherein force on the needle guard exceeding a predetermined amount overcomes a frictional force of the resistance band allowing the needle guard to retract relative to the body to unlock the grip relative to the body such that movement of the grip relative to the body causes medication to be delivered from the body.

**18 Claims, 36 Drawing Sheets**

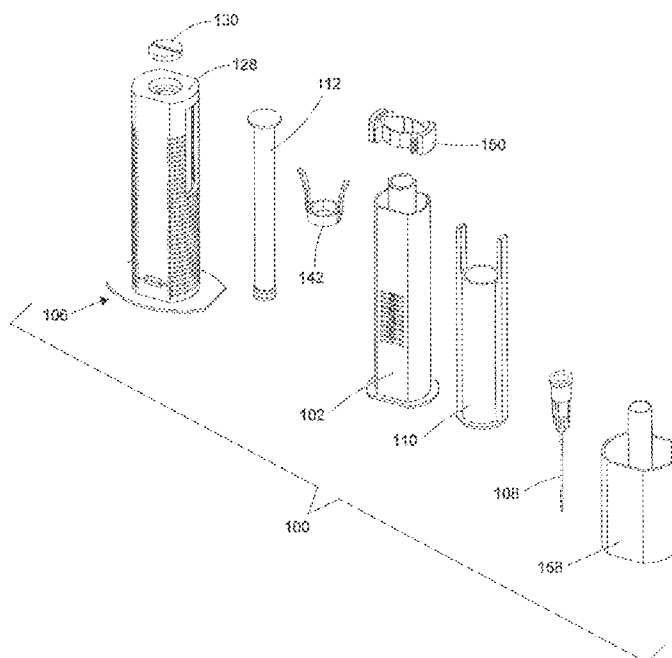

Related U.S. Application Data	(56)	References Cited
continuation-in-part of application No. 14/601,420, filed on Jan. 21, 2015, now abandoned.		U.S. PATENT DOCUMENTS
(60) Provisional application No. 61/929,768, filed on Jan. 21, 2014.	5,599,309 A *	2/1997 Marshall ..... A61M 5/46 604/117
	5,795,336 A *	8/1998 Romano ..... A61M 5/3271 604/110
(51) <b>Int. Cl.</b>	7,500,963 B2	3/2009 Westbye et al.
<i>A61M 5/31</i> (2006.01)	8,337,472 B2	12/2012 Edginton et al.
<i>A61M 5/32</i> (2006.01)	8,591,463 B1	11/2013 Cowe
(52) <b>U.S. Cl.</b>	2008/0228147 A1	9/2008 David-Hegerich et al.
CPC ..... <i>A61M 5/3287</i> (2013.01); <i>A61M 5/46</i>	2009/0259178 A1	10/2009 Brechbuchler et al.
(2013.01); <i>A61M 5/31501</i> (2013.01); <i>A61M</i>	2011/0092915 A1 *	4/2011 Olson ..... A61M 5/2033 604/198
<i>2005/3114</i> (2013.01); <i>A61M 2005/3126</i>	2011/0202011 A1	8/2011 Wozencroft
(2013.01); <i>A61M 2005/3139</i> (2013.01); <i>A61M</i>	2012/0203186 A1	8/2012 Vogt et al.
<i>2005/3143</i> (2013.01); <i>A61M 2205/586</i>	2012/0233834 A1	9/2012 Szechinski et al.
(2013.01)	2012/0316508 A1 *	12/2012 Kirchhofer ..... A61M 5/31553 604/198
(58) <b>Field of Classification Search</b>	2013/0204229 A1	8/2013 Olson et al.
CPC ..... A61M 5/31511; A61M 5/31501; A61M	2013/0237921 A1	9/2013 Lannan et al.
2005/3126; A61M 2005/3139	2013/0317480 A1	11/2013 Reber et al.
See application file for complete search history.	2017/0165429 A1	6/2017 Holmqvist

\* cited by examiner

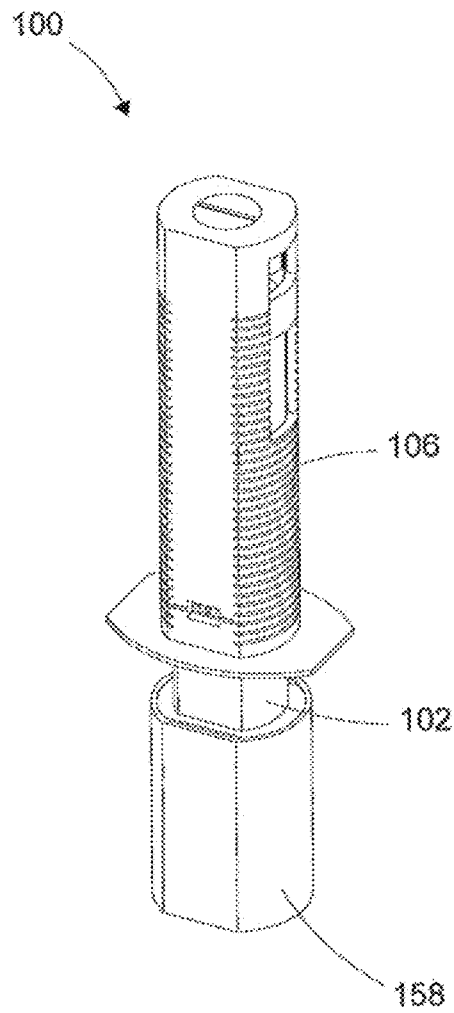
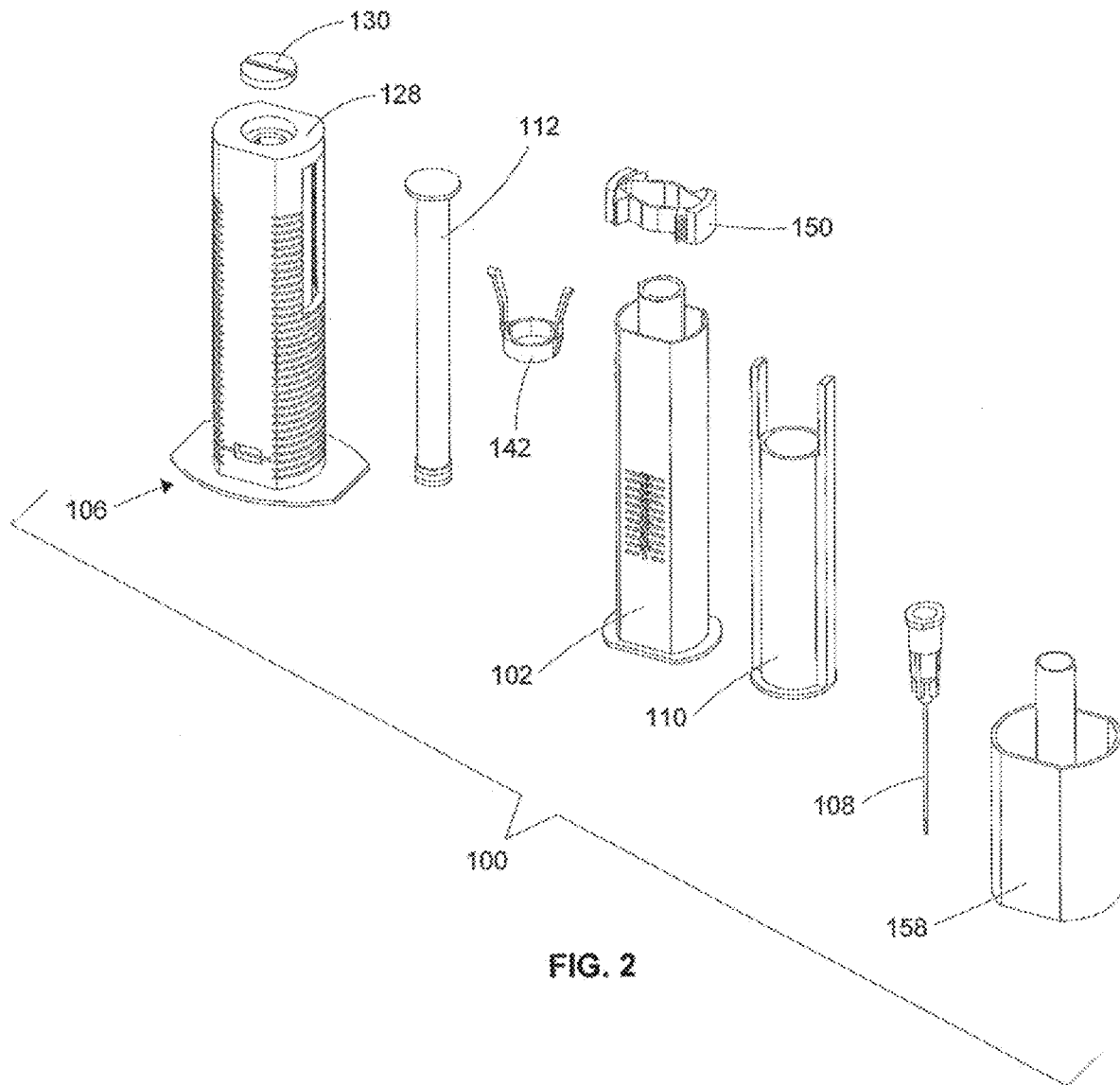


FIG. 1



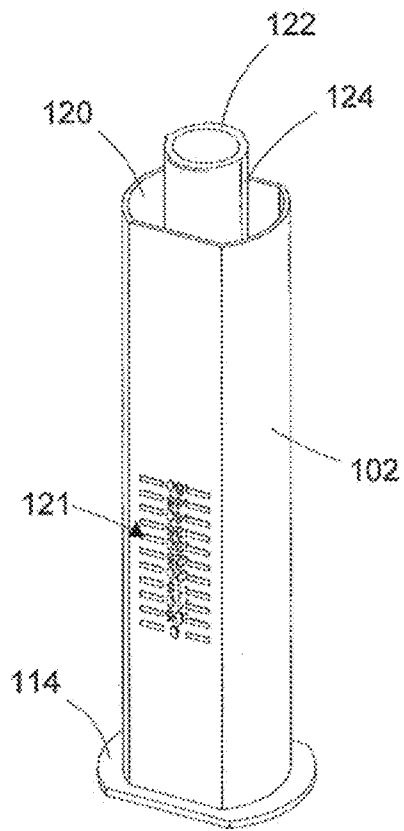


FIG. 3

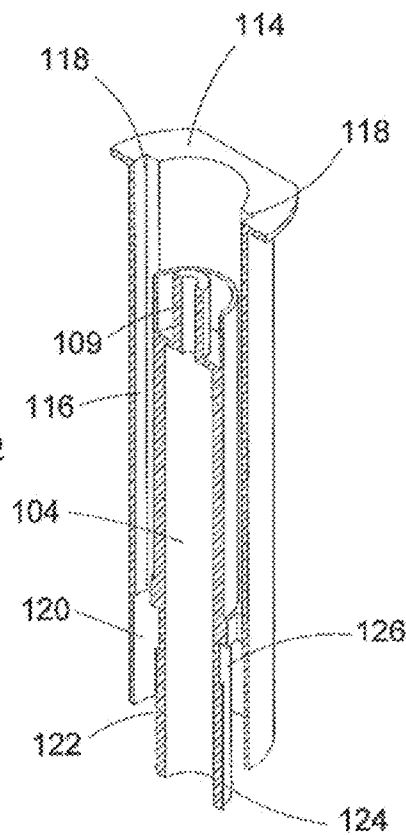


FIG. 4

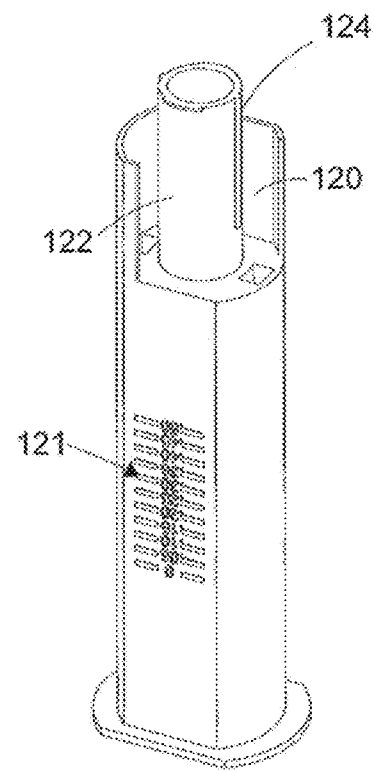


FIG. 5

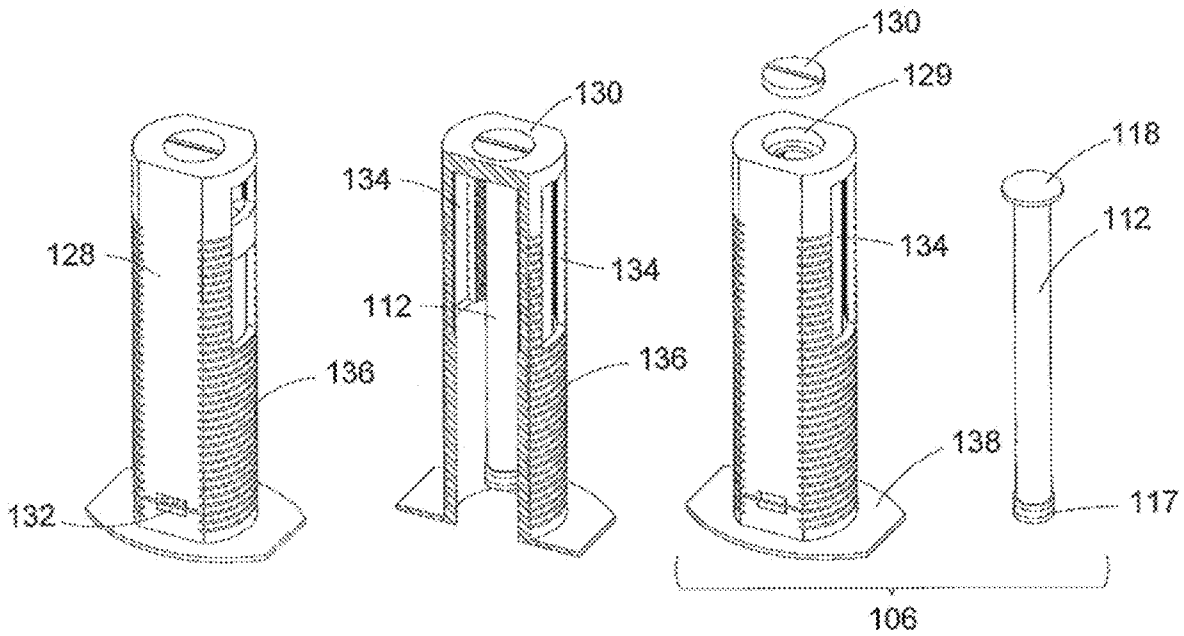


FIG. 6

FIG. 7

FIG. 8

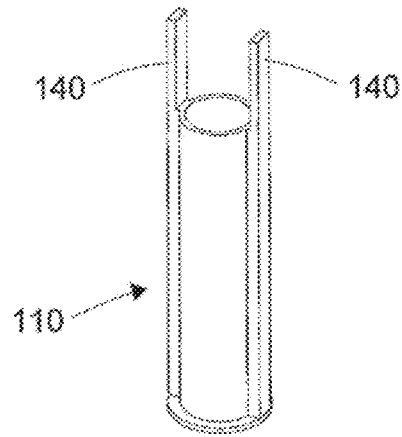


FIG. 9

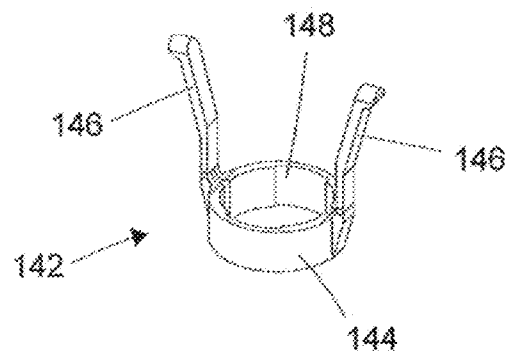


FIG. 10

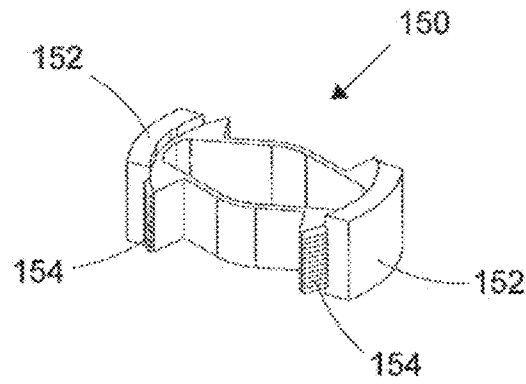


FIG. 11

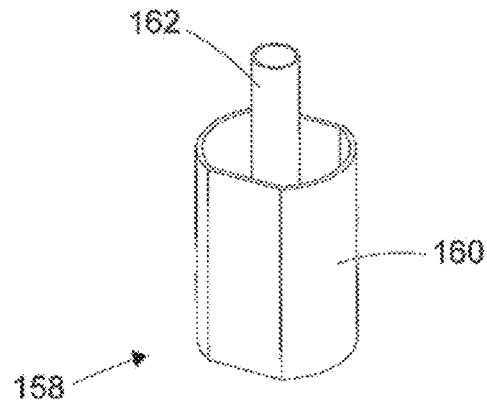


FIG. 12



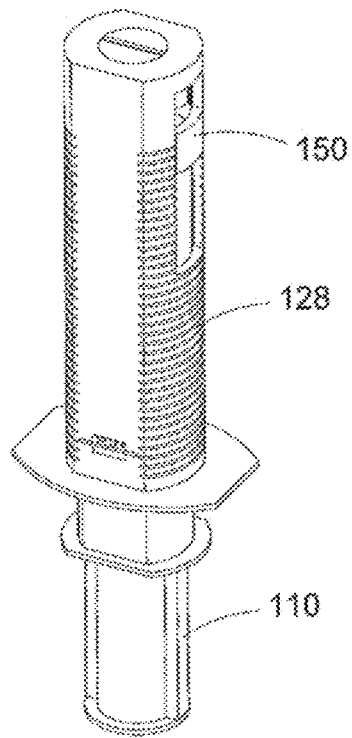


FIG. 13

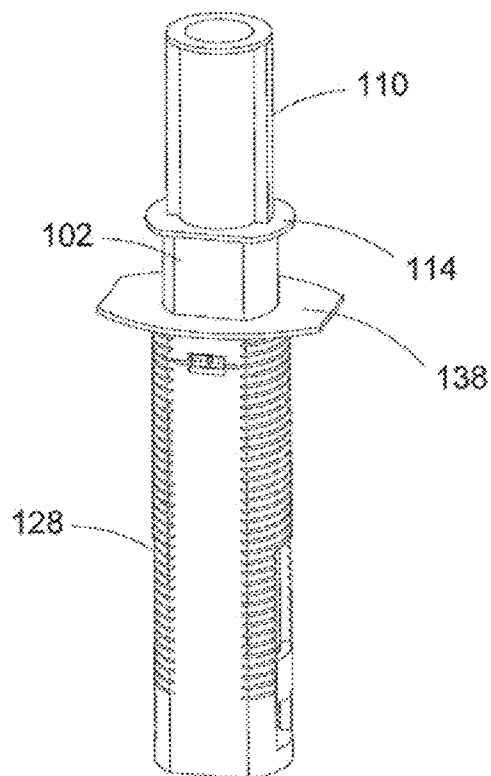


FIG. 14

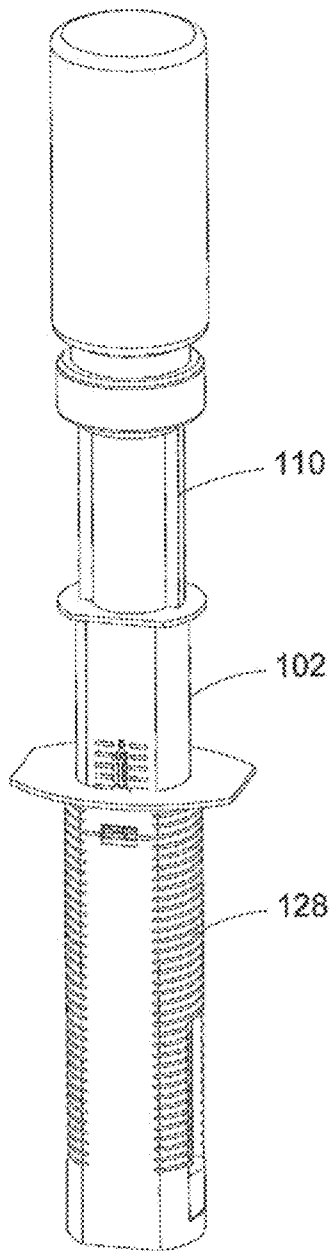


FIG. 15

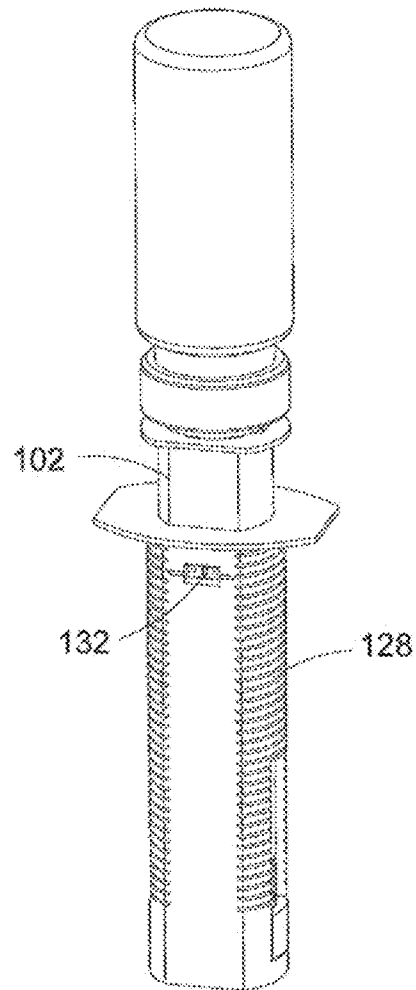


FIG. 16

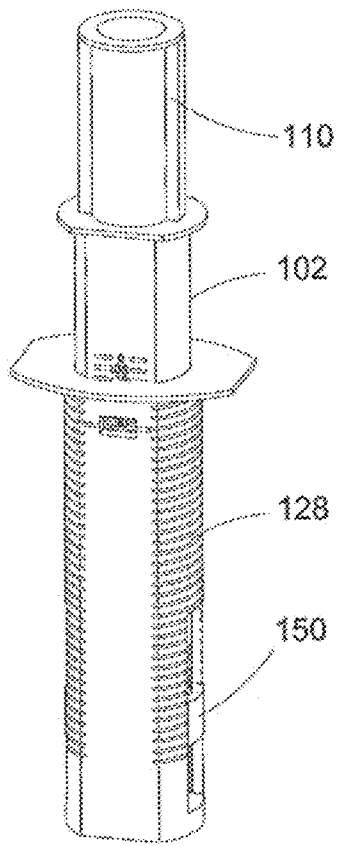


FIG. 17

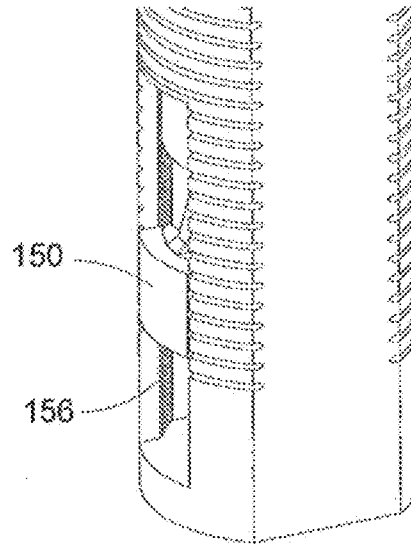


FIG. 18

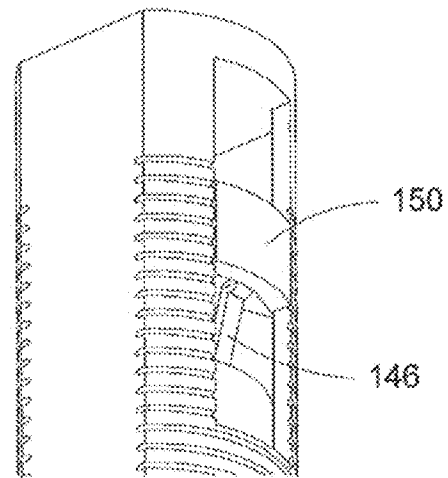


FIG. 19

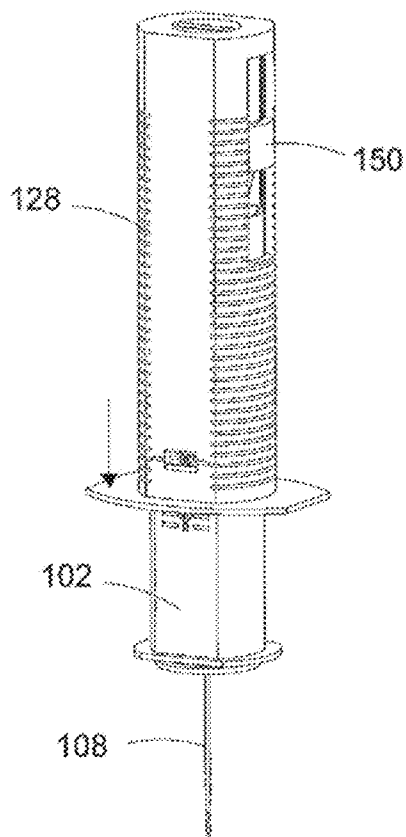


FIG. 20

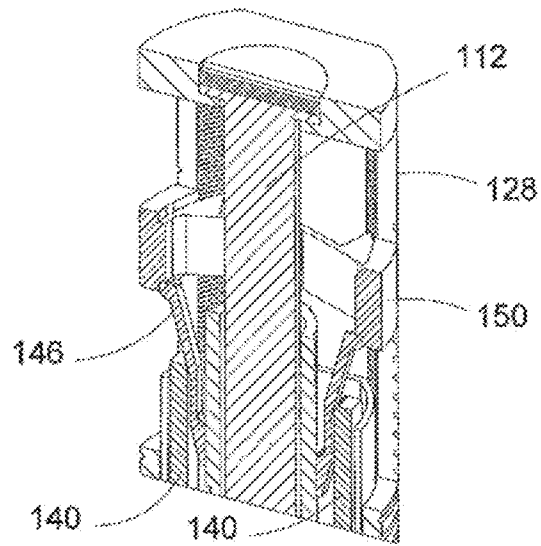


FIG. 21

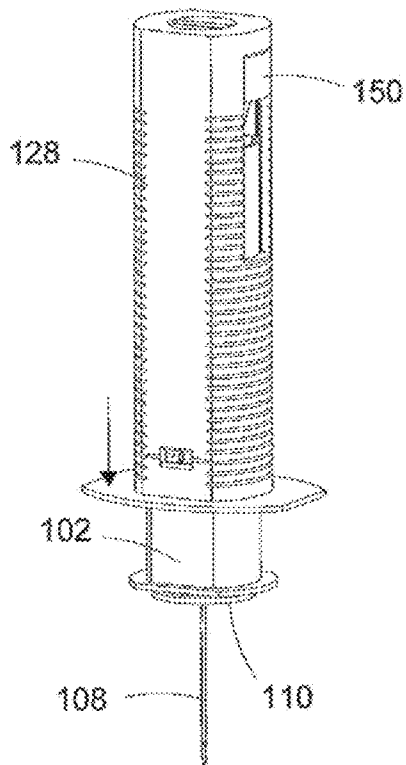


FIG. 22

200

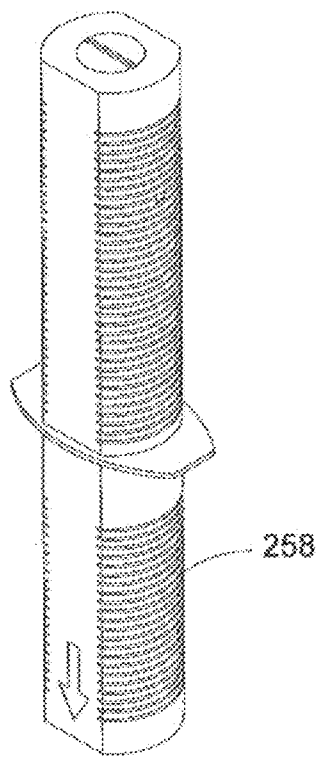


FIG. 23

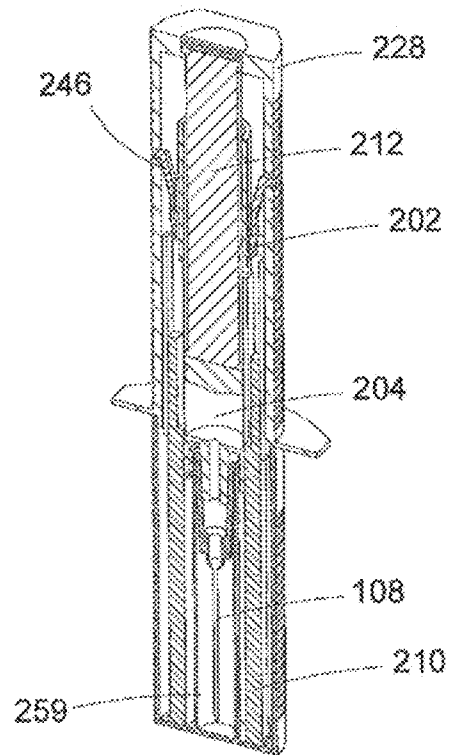


FIG. 24

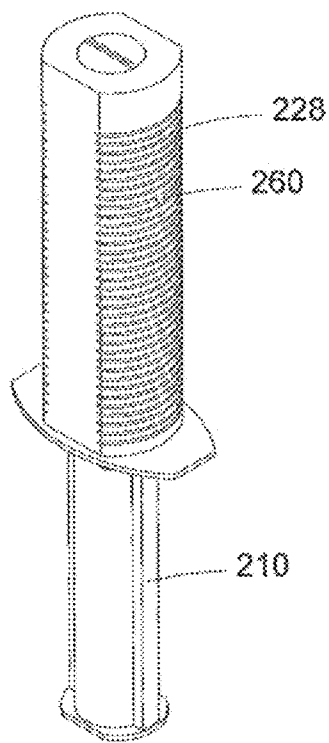


FIG. 25

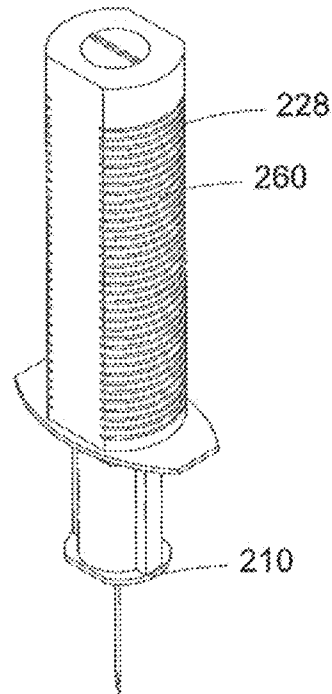
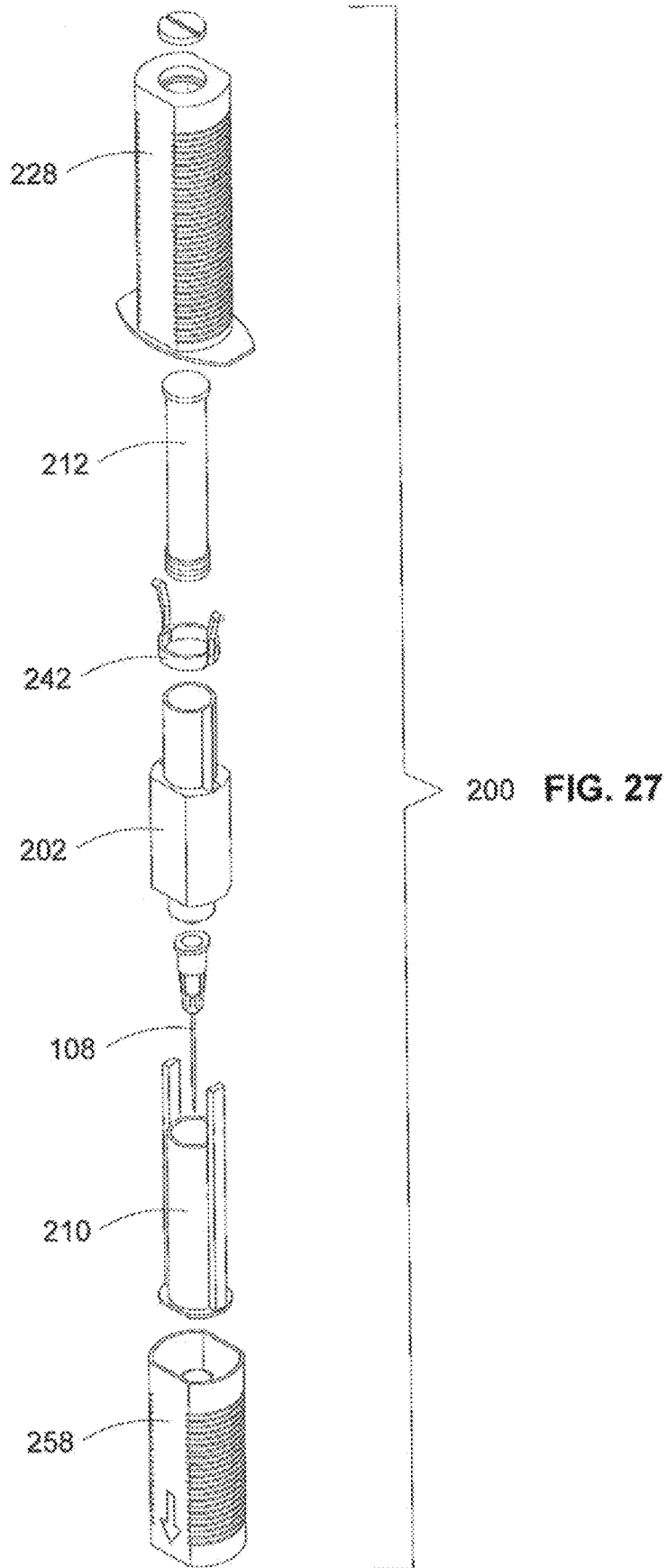
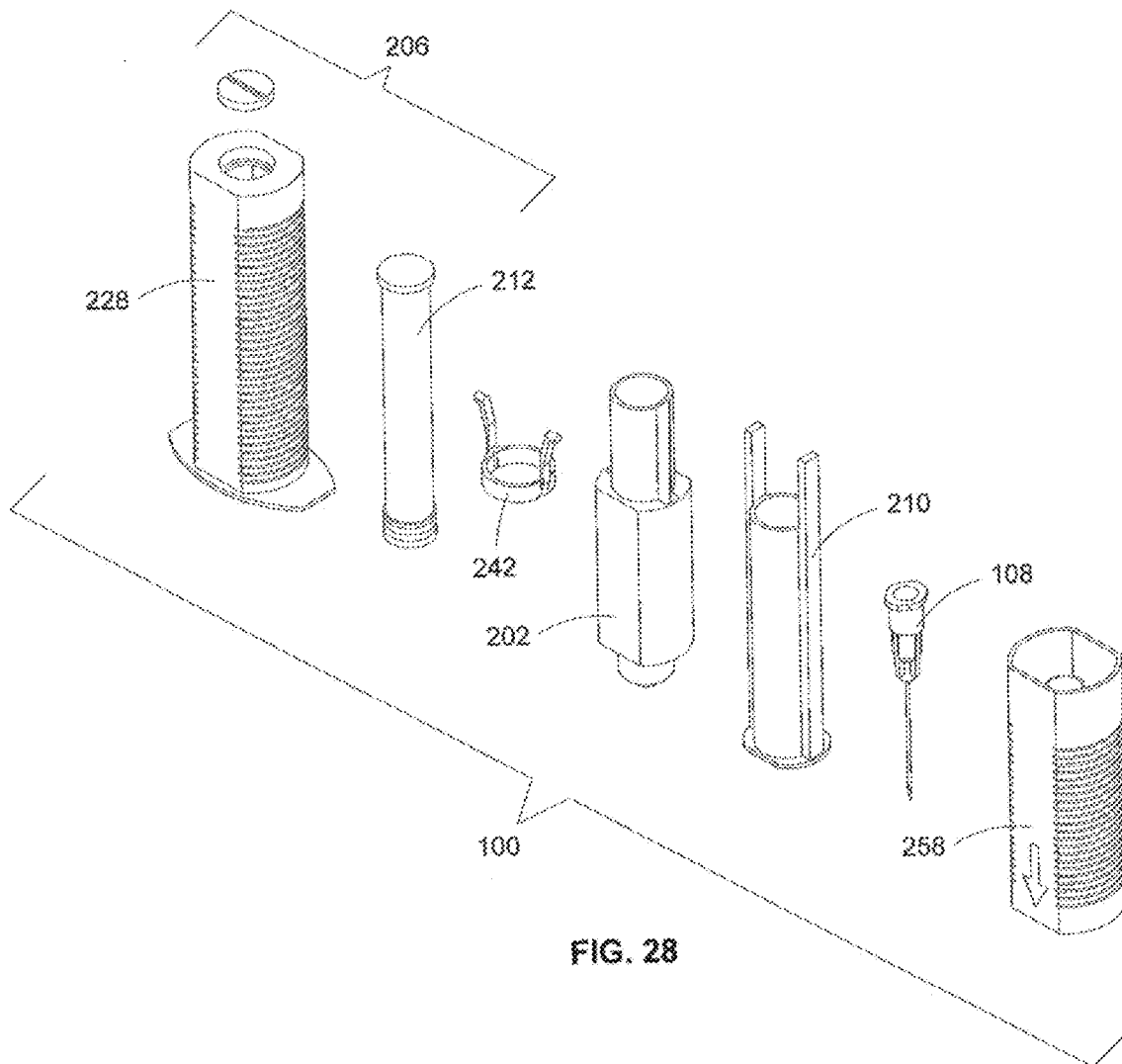


FIG. 26







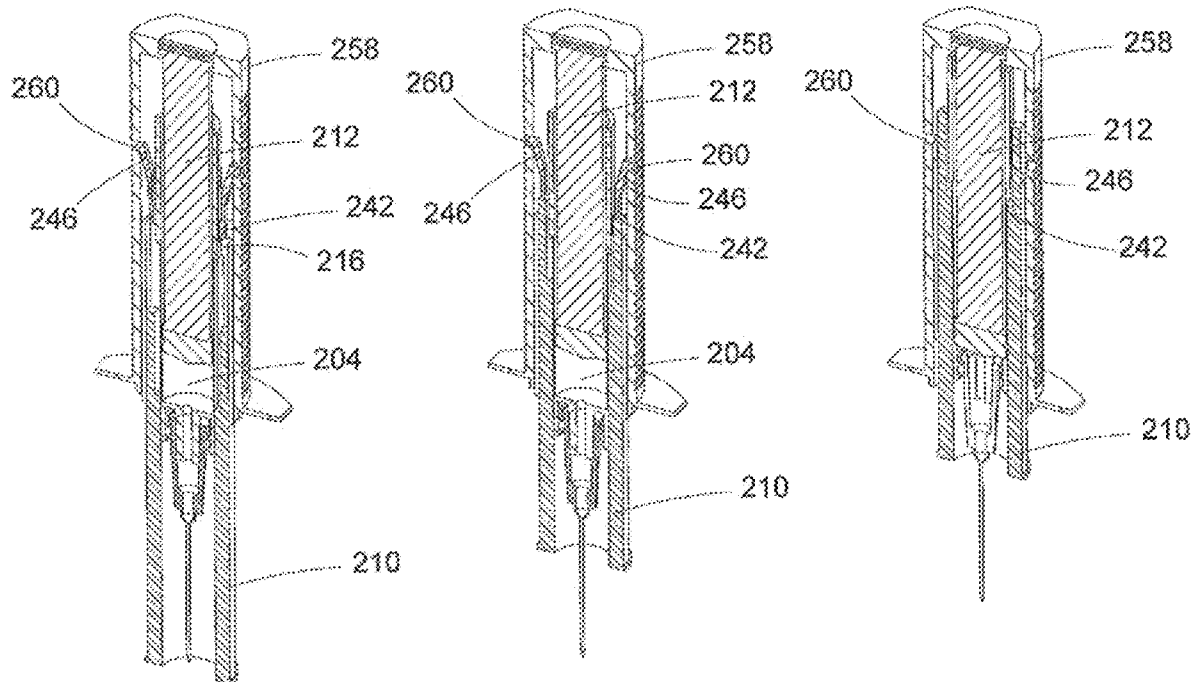


FIG. 29

FIG. 30

FIG. 31

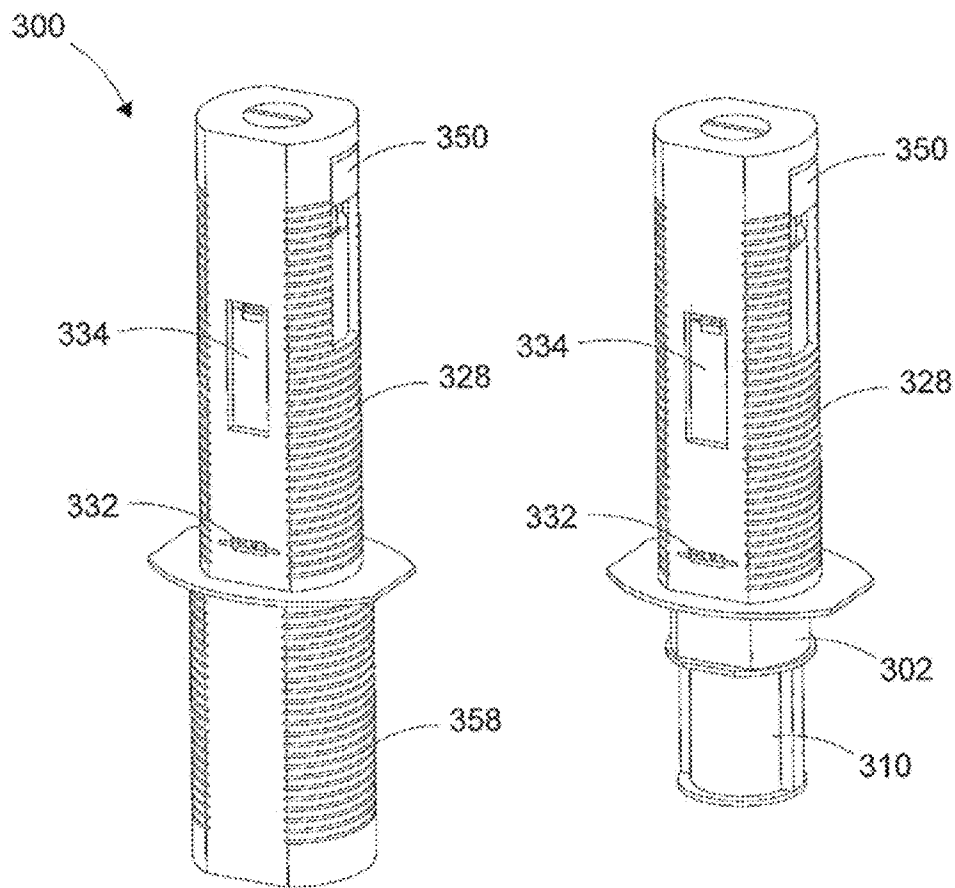


FIG. 32

FIG. 33

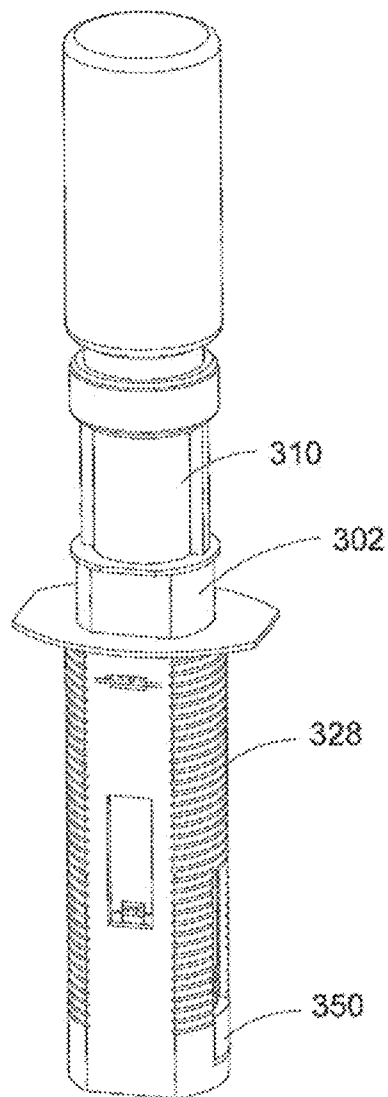


FIG. 34

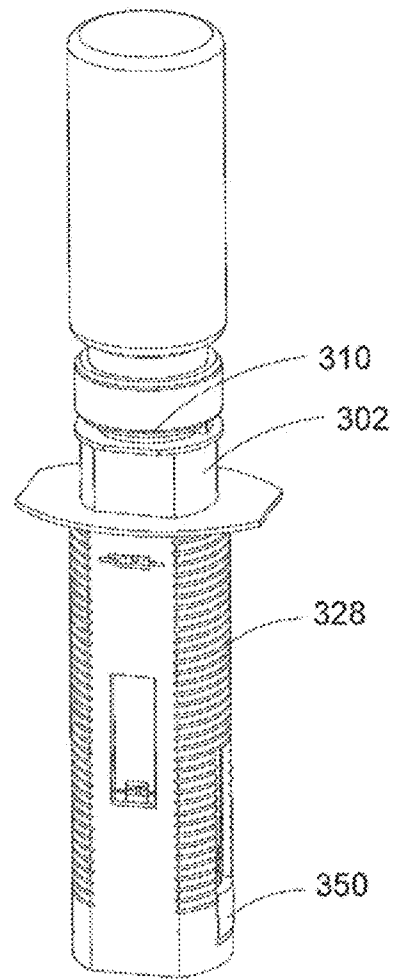


FIG. 35

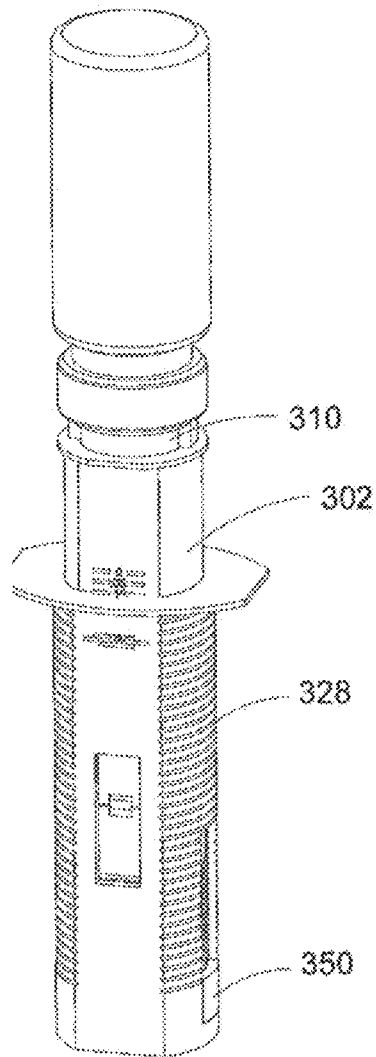


FIG. 36

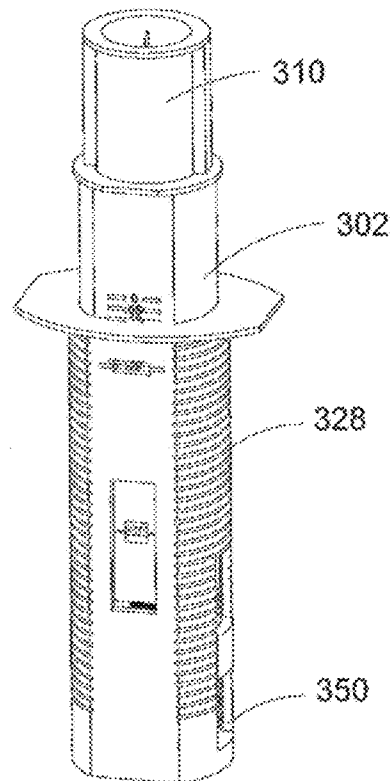


FIG. 37

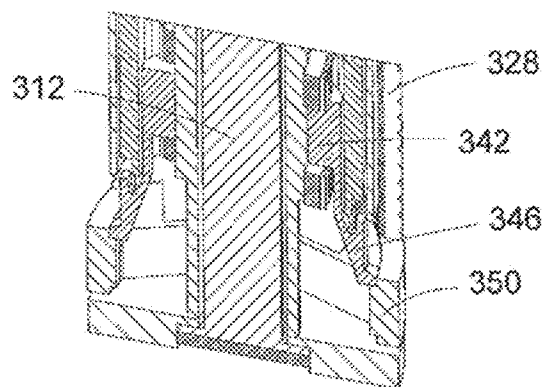


FIG. 38

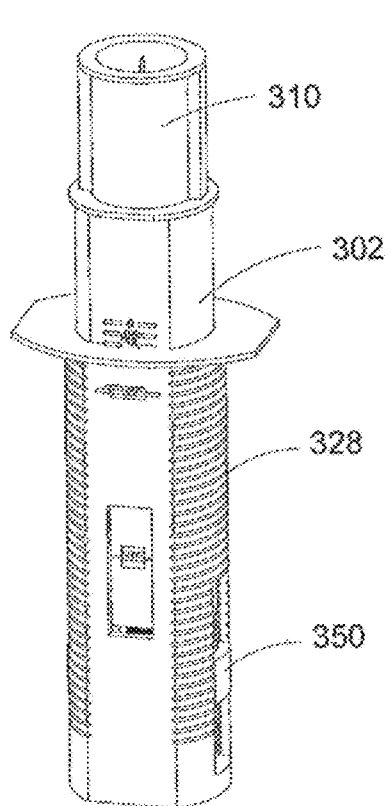


FIG. 39

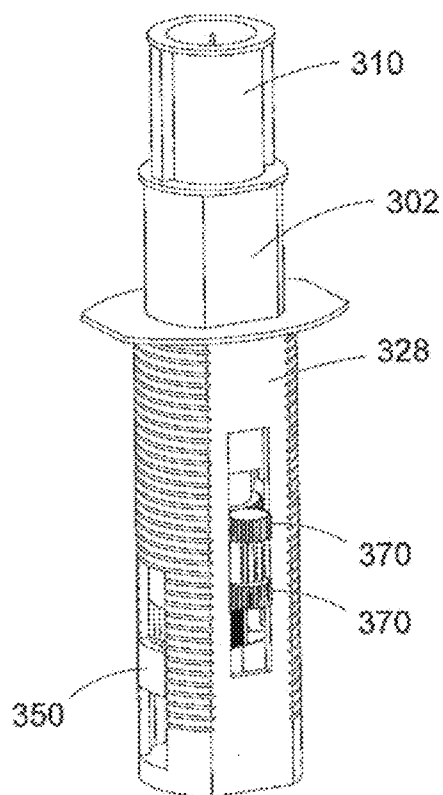


FIG. 40

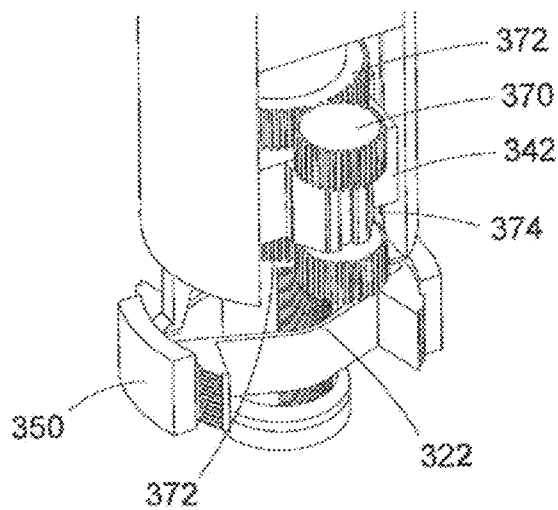


FIG. 41

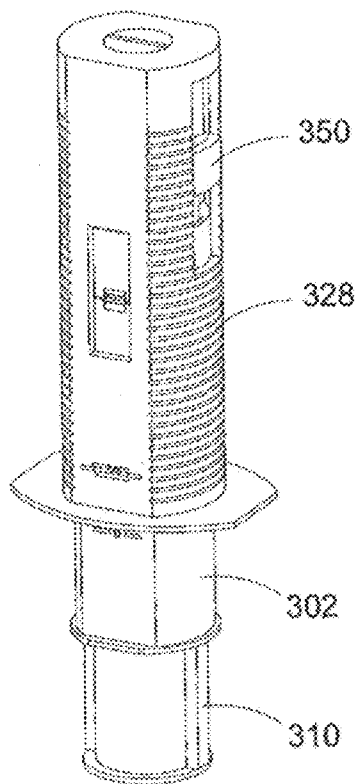


FIG. 42

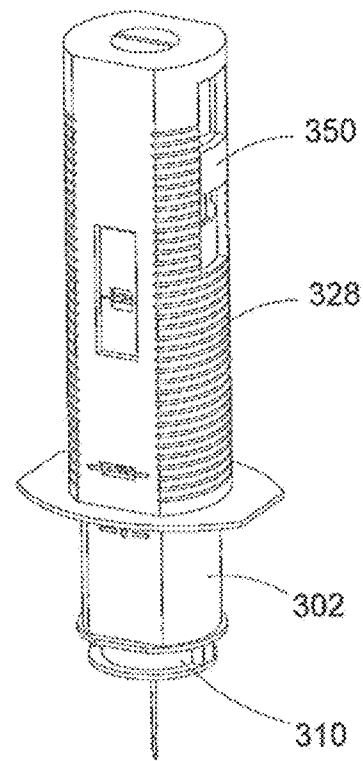


FIG. 43

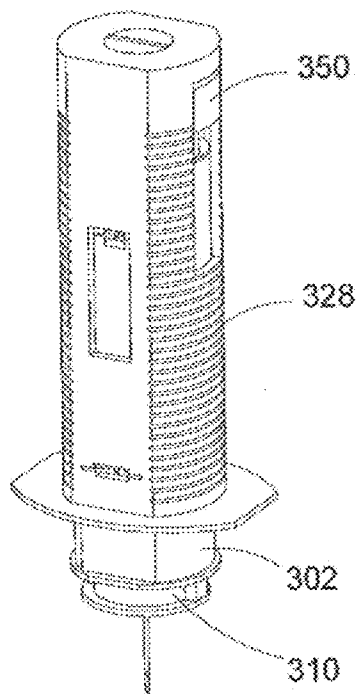


FIG. 44

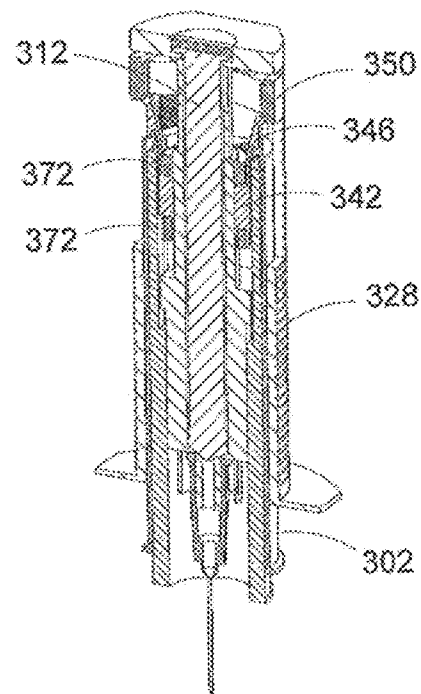


FIG. 45

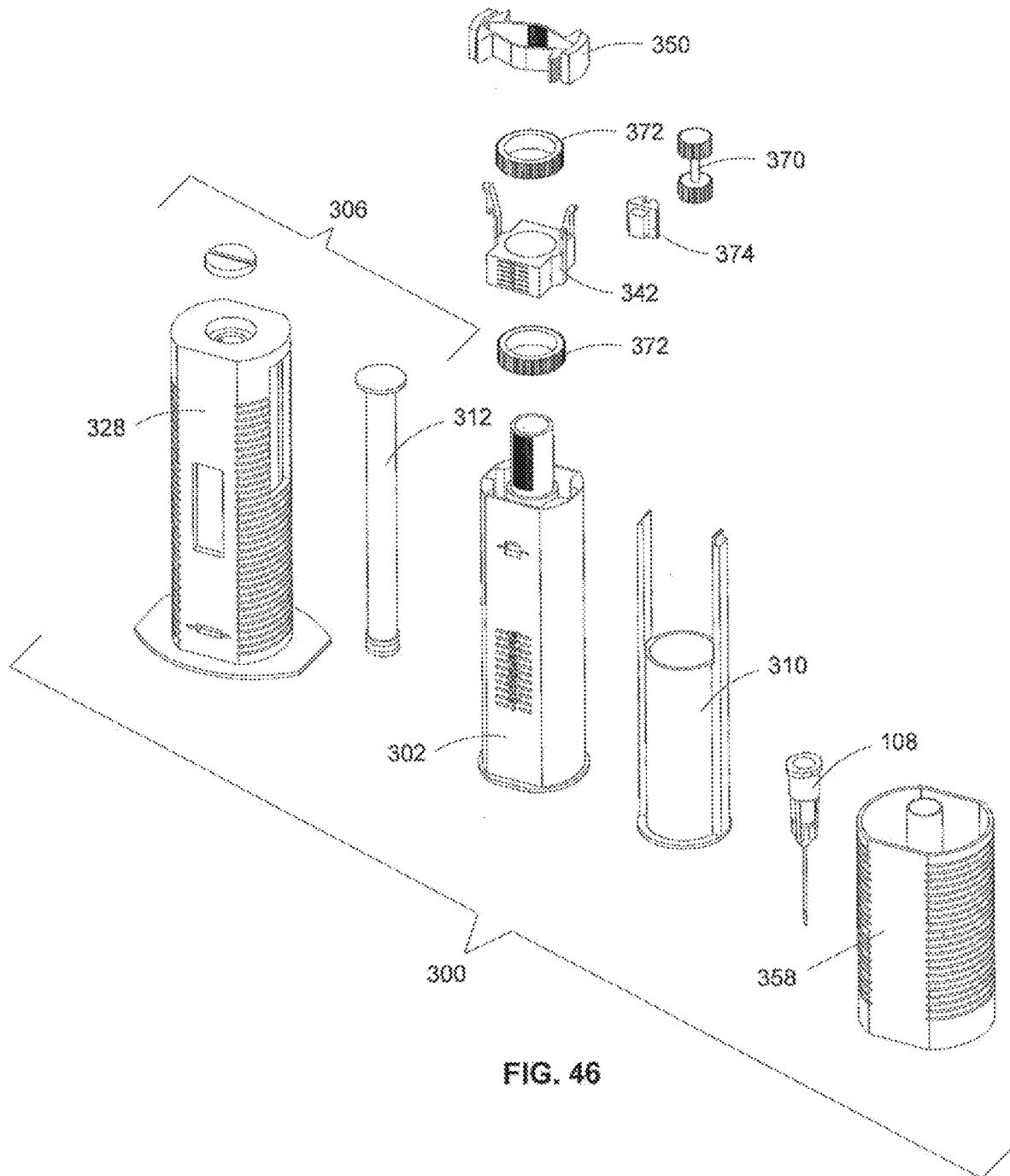


FIG. 46

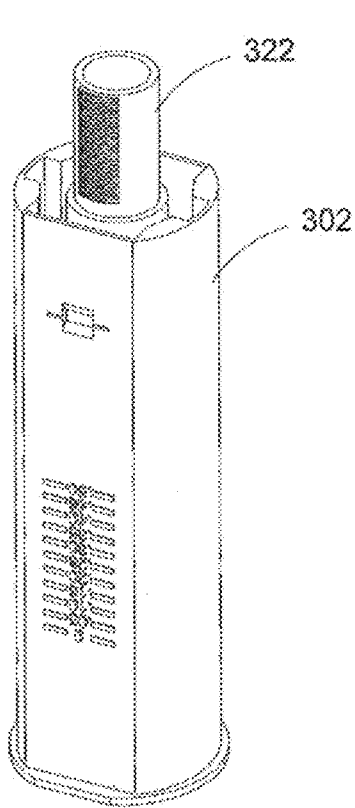


FIG. 47

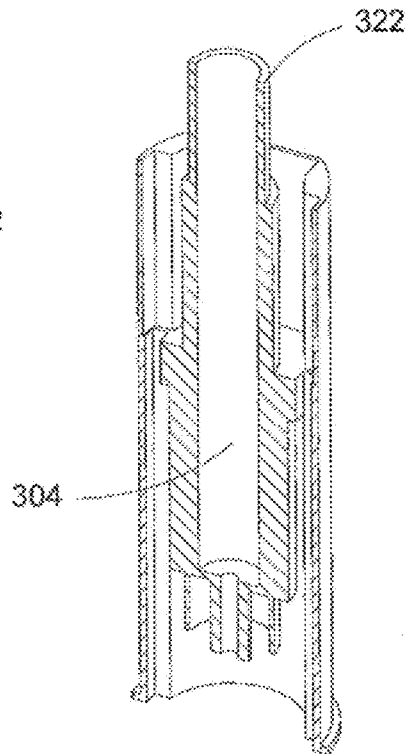


FIG. 48

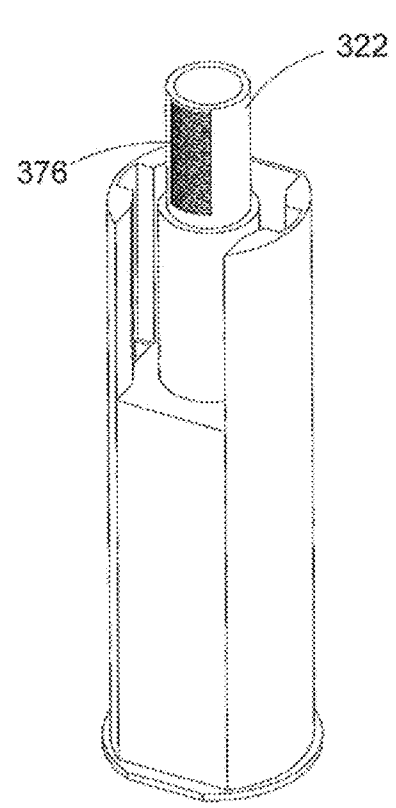


FIG. 49



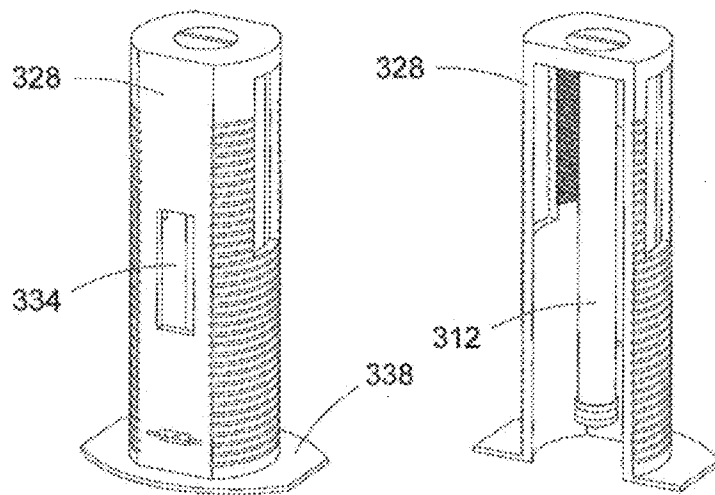


FIG. 50

FIG. 51

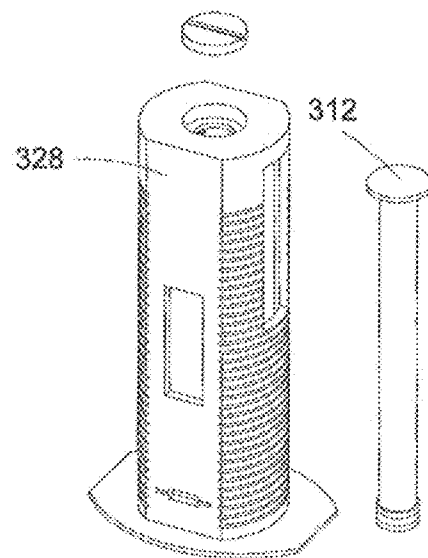


FIG. 52

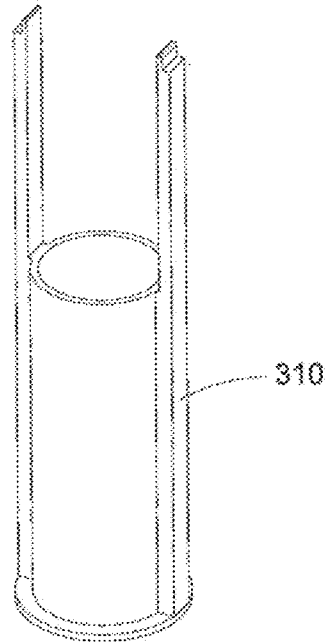


FIG. 53

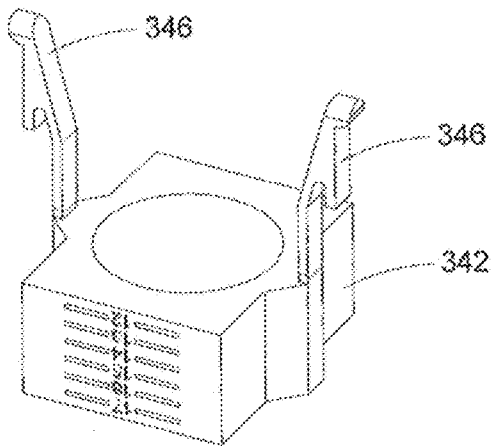


FIG. 54

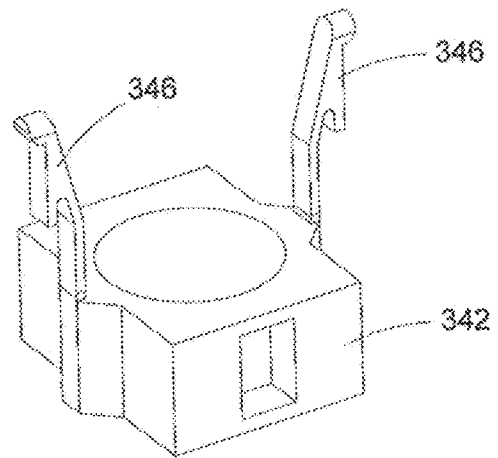


FIG. 55

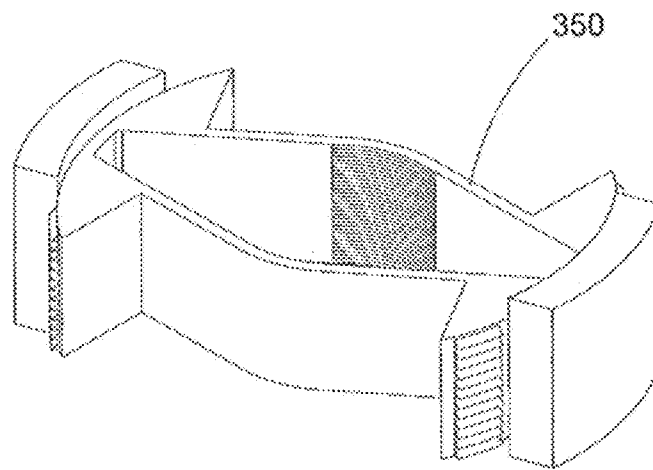


FIG. 56

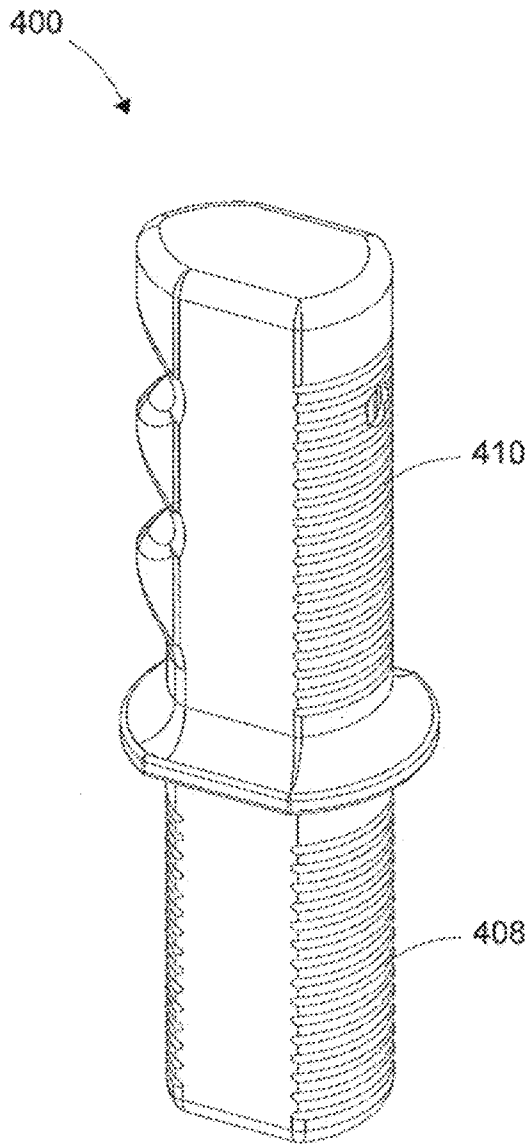


FIG. 57

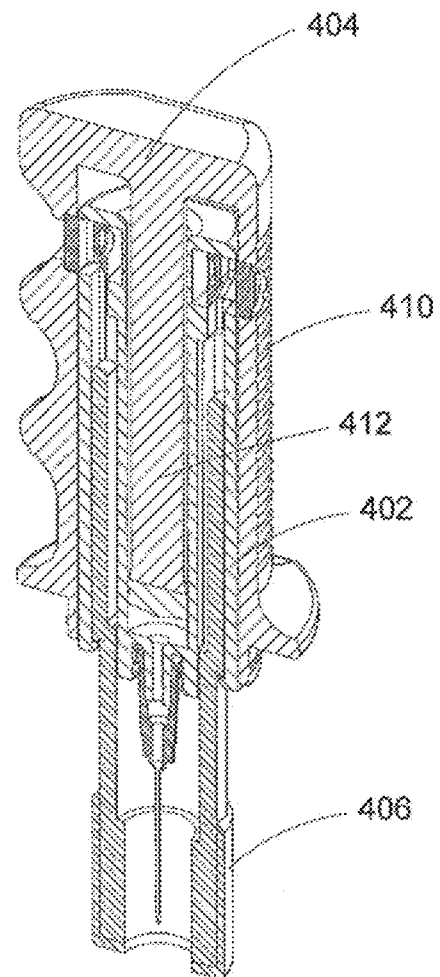


FIG. 58

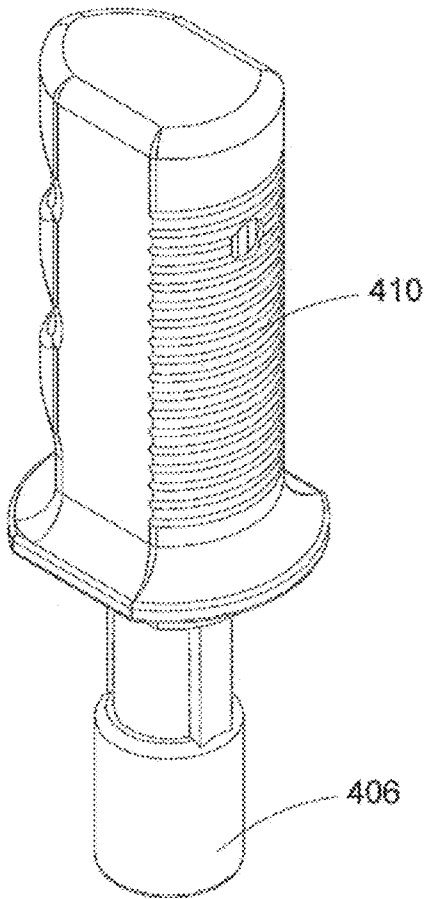


FIG. 59

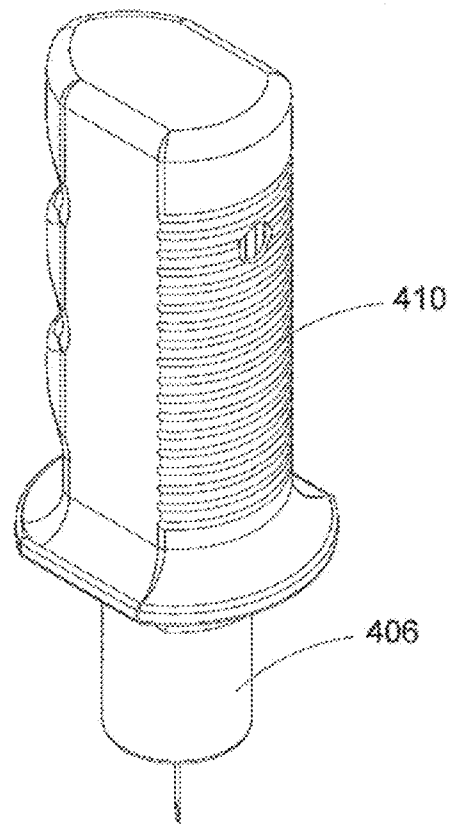


FIG. 60

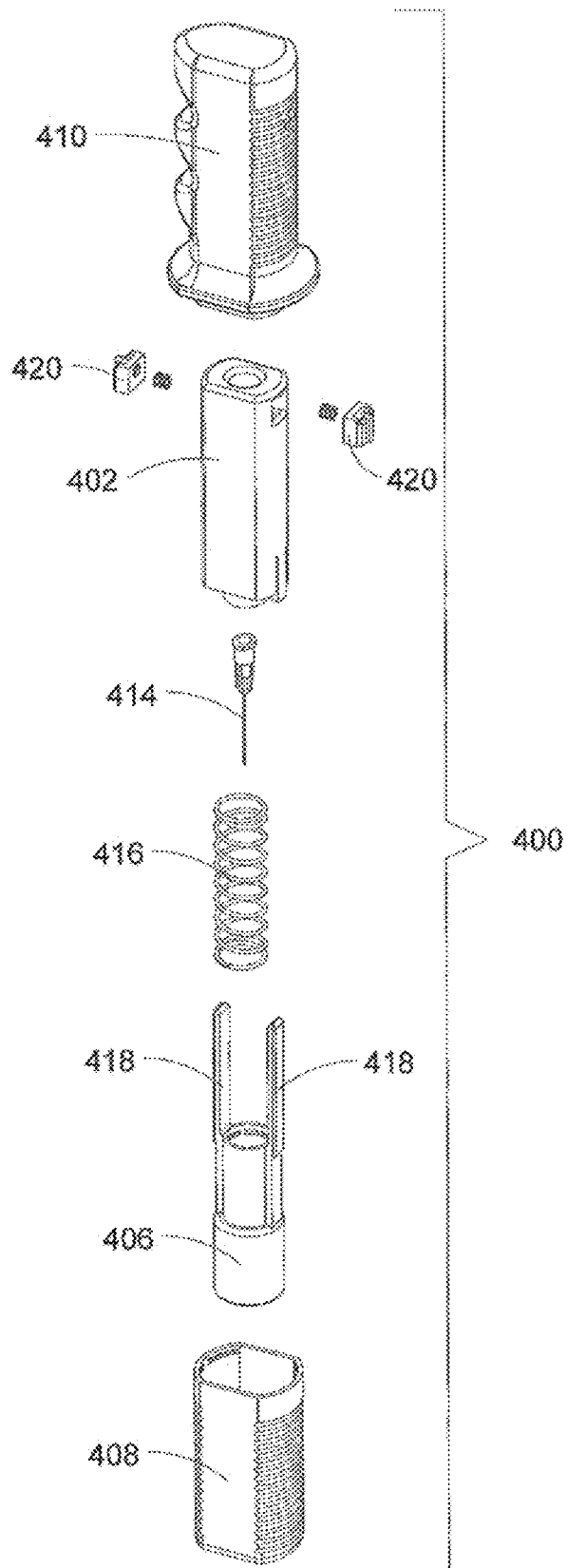
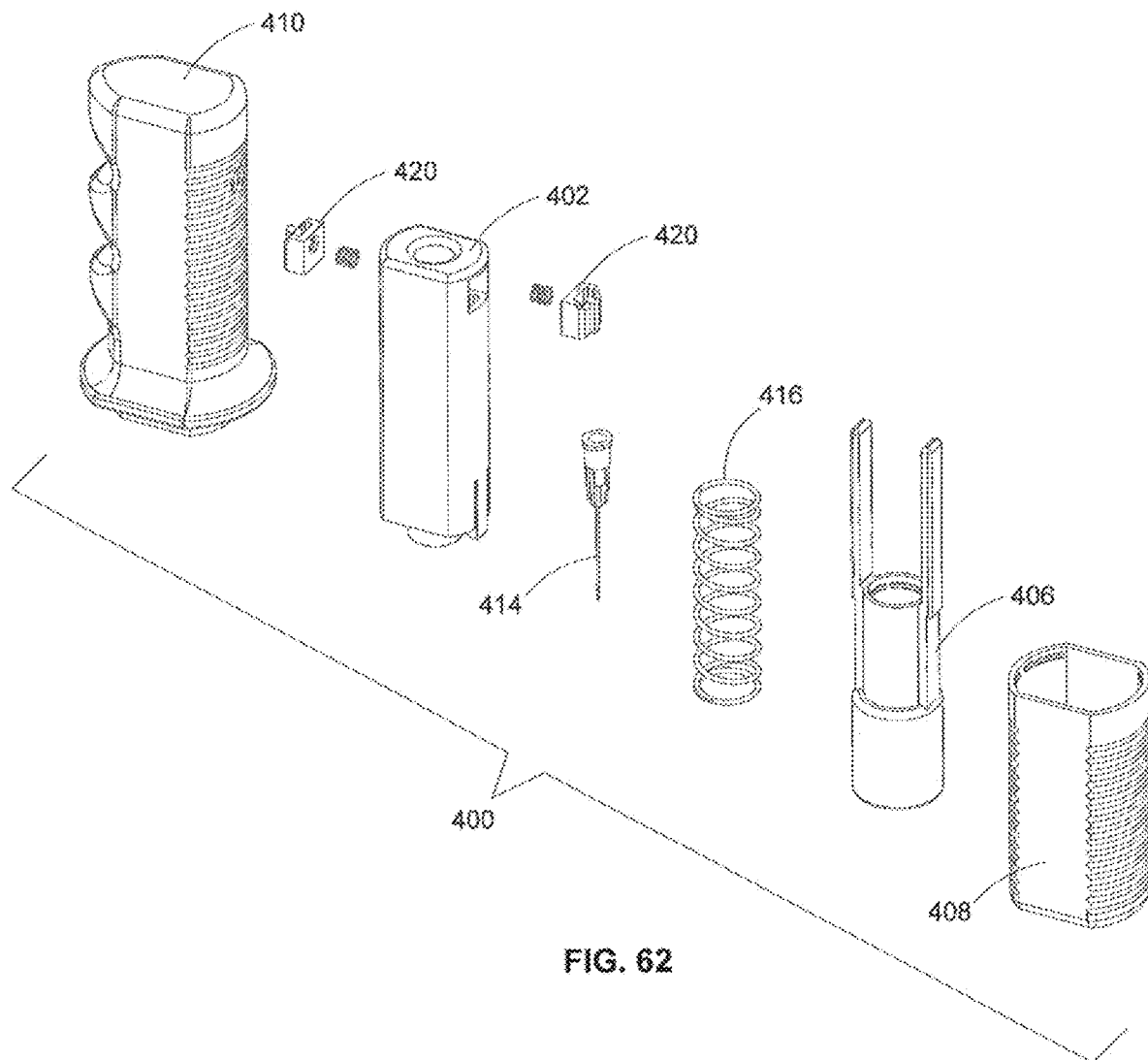


FIG. 61



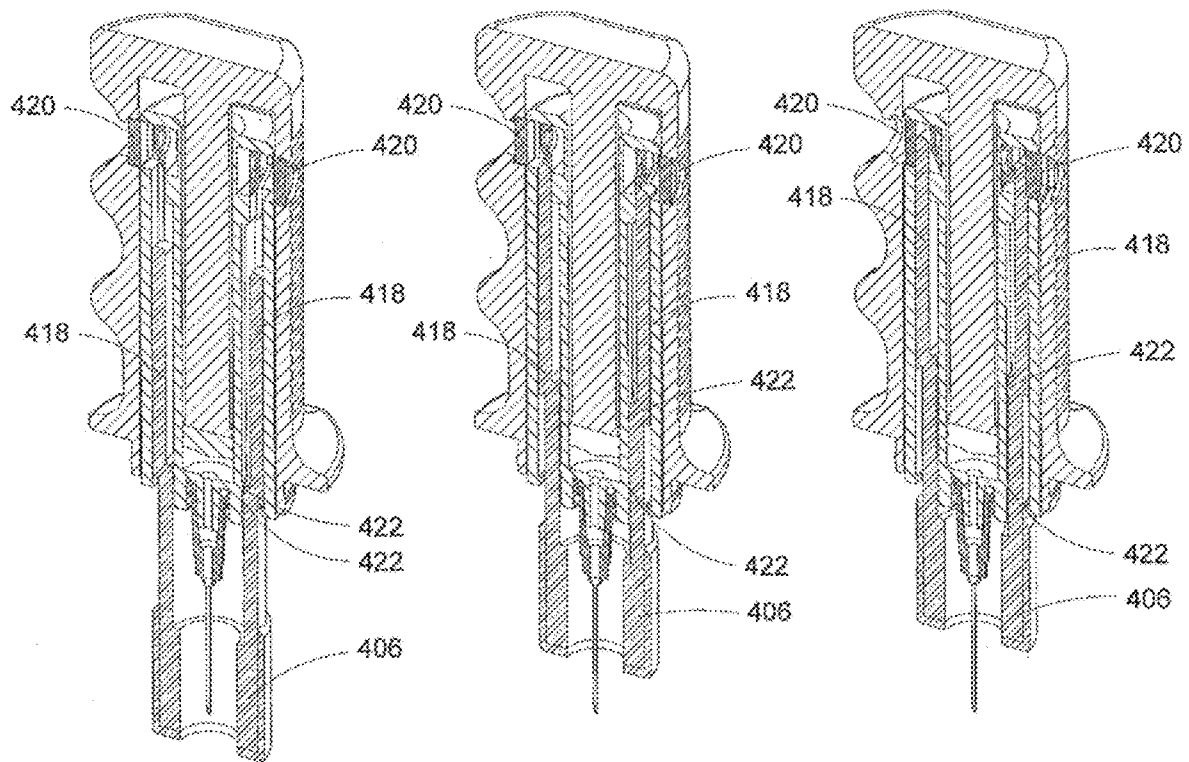


FIG. 63

FIG. 64

FIG. 65



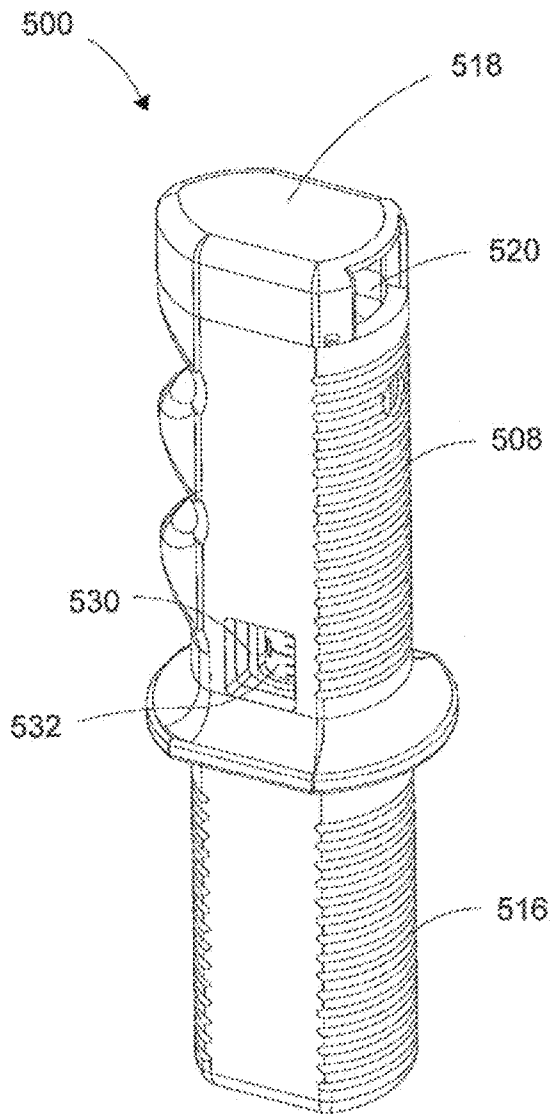


FIG. 66

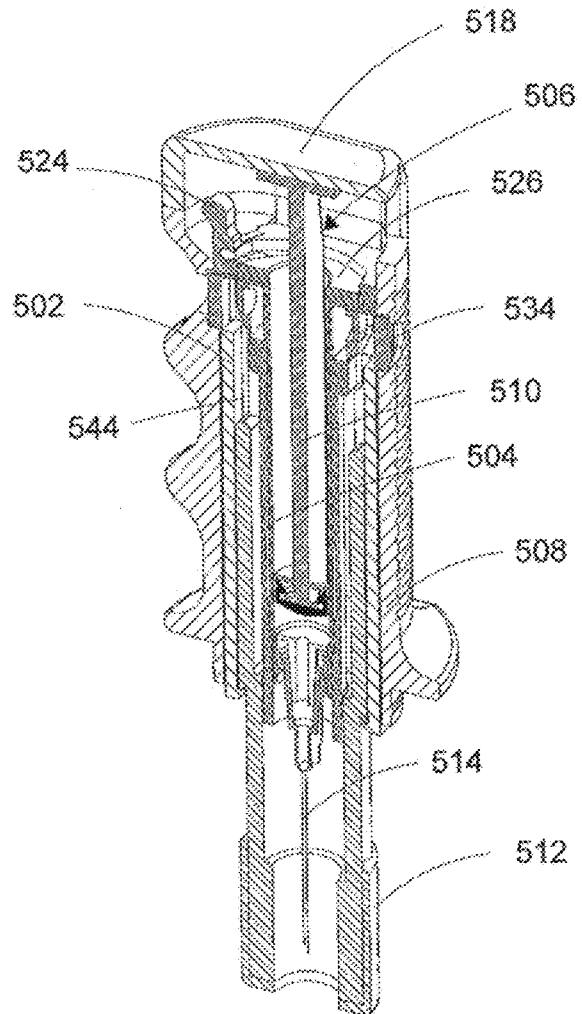


FIG. 67

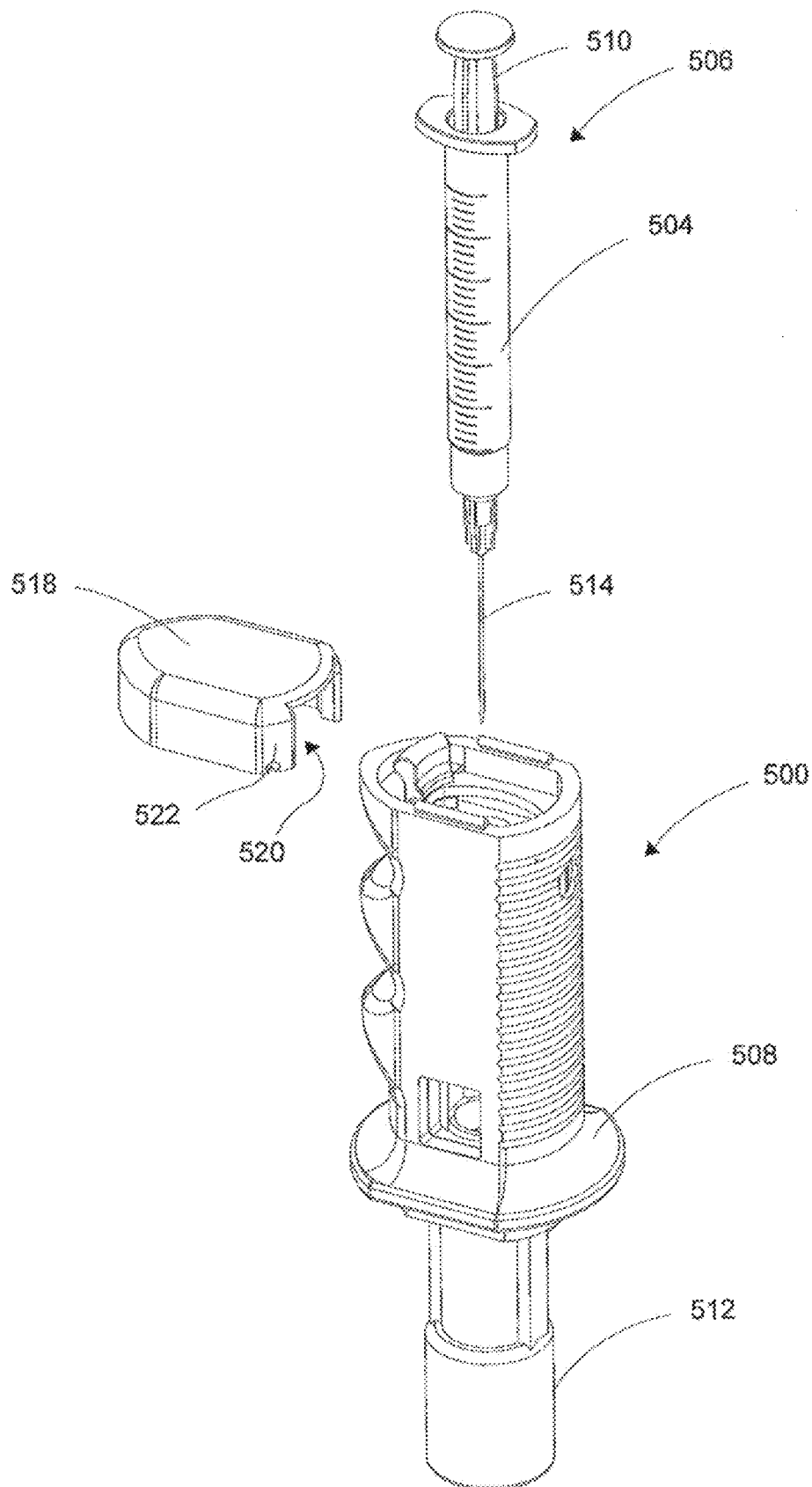


FIG. 68

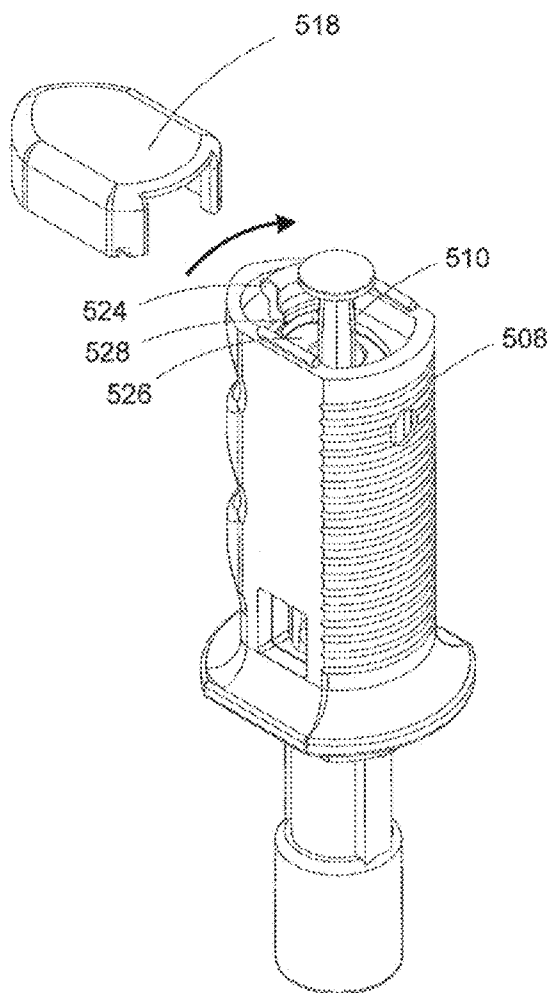


FIG. 69

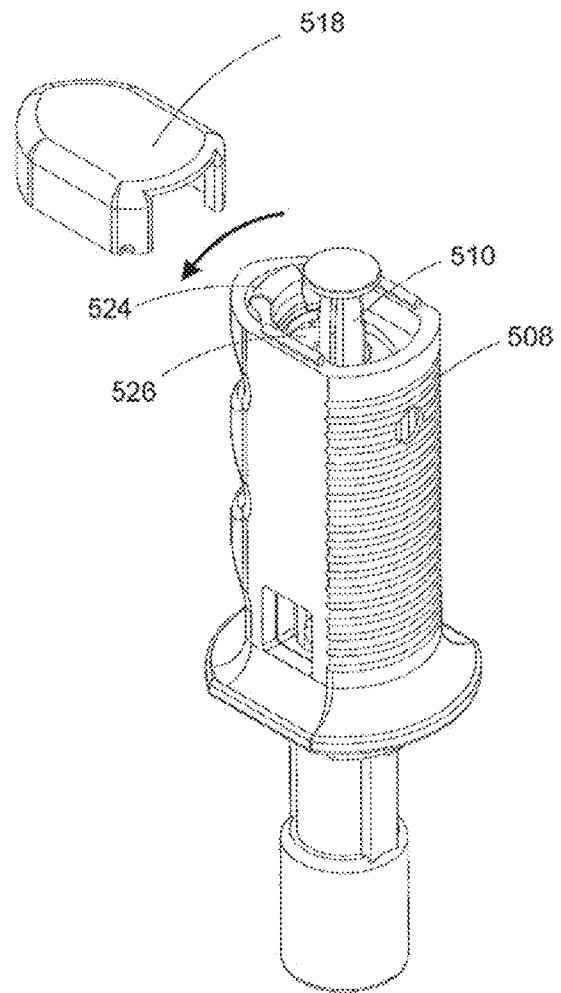


FIG. 70

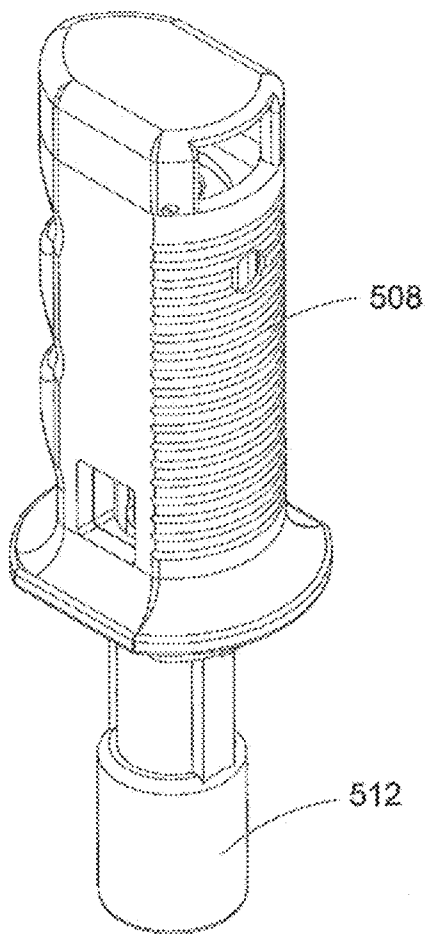


FIG. 71

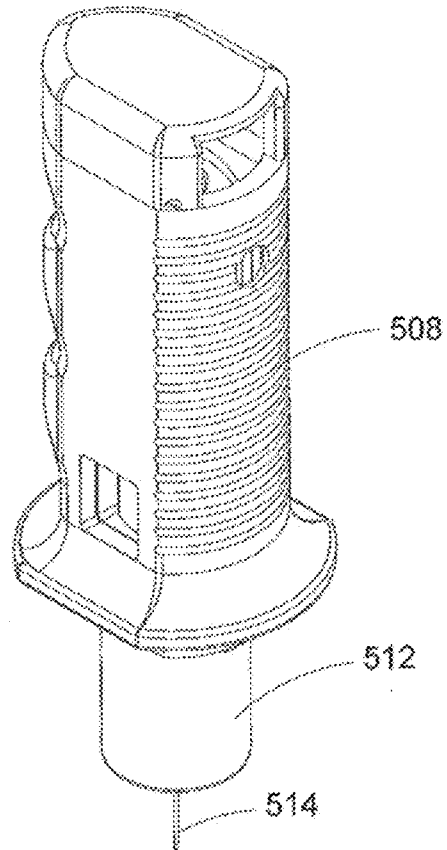


FIG. 72


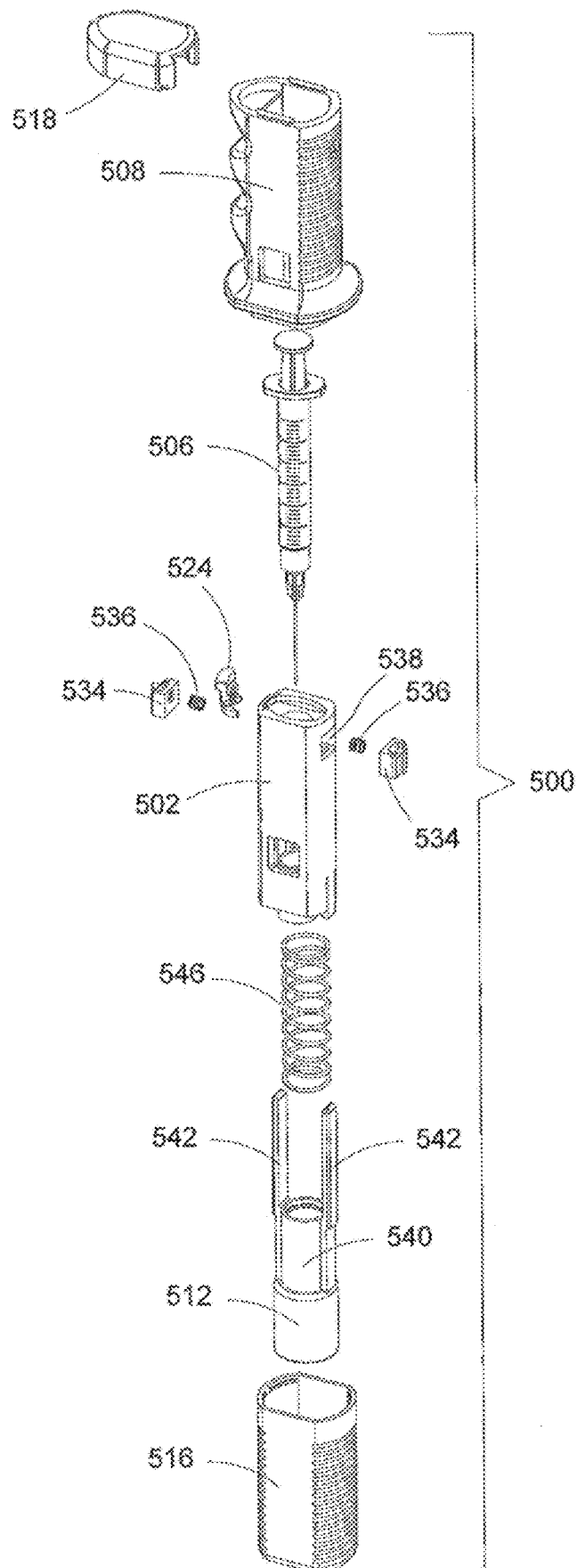


FIG. 73

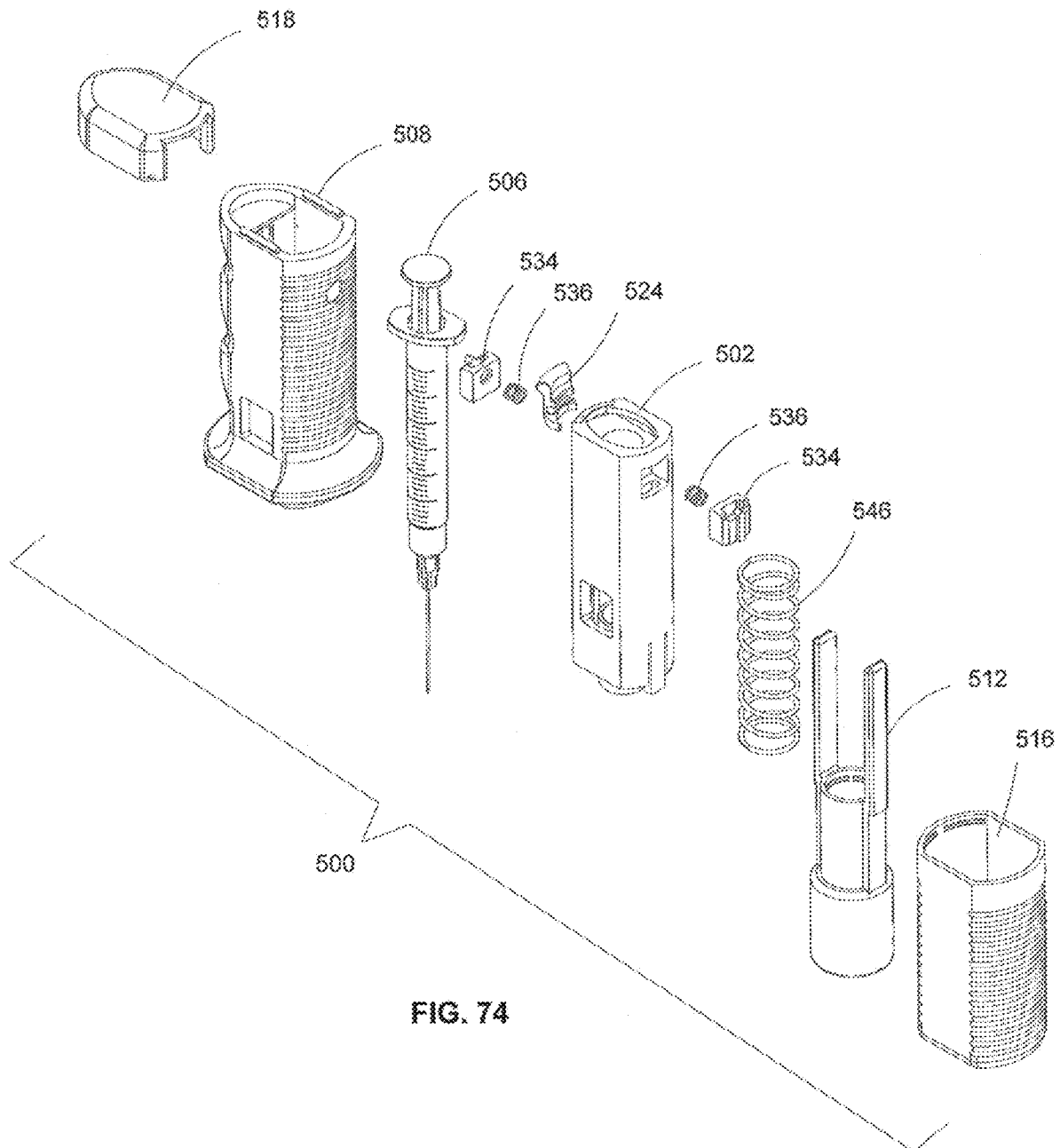

1

**FORCE ACTUATED INJECTION DEVICE****CROSS-REFERENCE TO RELATED APPLICATIONS**

This Continuation application claims priority to U.S. application Ser. No. 15/200,298 filed Jul. 1, 2016, which is a Continuation-in-part of U.S. application Ser. No. 14/601,420 filed Jan. 21, 2015 claiming priority from U.S. App. No. 61/929,768 filed Jan. 21, 2014, each of which is incorporated by reference herein in its entirety.

**TECHNICAL FIELD AND BACKGROUND OF THE INVENTION**

The present invention relates generally to an injection device for improving the speed, accuracy, safety and effectiveness of the administration of an emergency, life-saving injection during an anaphylactic event, and more particularly, to an injection device that operates by way of sequential manual actuation, wherein movement of a needle guard in response to pressure against the injection site initiates a hypodermic injection and triggers the release of a plunger for delivering medication from a reservoir, in one continuous motion.

Anaphylactic reactions are serious and sometimes fatal reactions to allergens, most commonly caused by food, drugs, insect stings, etc. Studies suggest that food allergies are increasing world-wide. According to a study released by the CDC, food allergies increased approximately 50% between 1997 and 2011. More than 17 million Europeans have a food allergy according to the European Academy of Allergy and Clinical Immunology. The CDC has also reported that food allergies result in more than 300,000 ambulatory care visits per year among children under the age of 18.

Auto-injectors are universally recommended as first-aid treatment for anaphylaxis. In a World Allergy Organization survey conducted in 2003, auto-injectors containing 0.3 mg of epinephrine were reported to be available in 56.4% of countries, and those containing 0.15 mg of epinephrine were reported to be available in 43.6% of countries. Auto-injectors containing an infant dose were not available in any country.

There have been many studies published about the inability of auto-injectors to effectively deliver epinephrine, particularly in pediatric patients. A first disadvantage of auto injectors is the ineffectiveness to deliver the accurate amount of epinephrine during an emergency anaphylaxis. Current epinephrine auto injectors are available in two standard doses: 0.15 mg and 0.3 mg, with the 0.15 mg dose for patients 30 kg and under, and the 0.30 mg dose for patients above 30 kg. According to The Journal of Allergy and Clinical Immunology, neither of these doses is appropriate for children weighing less than 10 Kg. A 20 kg child would have a similar issue in being under dosed with the 0.15 mg dose, and overdosed with the 0.30 mg dose. An overdose or under dose during a life threatening anaphylaxis, especially for a child or patient with small body mass, could be fatal.

A second disadvantage of auto injectors is the inability to effectively deliver the medication intramuscularly using one standard needle length. A study conducted at the Phoenix Children's Hospital found that the needle on epinephrine auto-injectors is not long enough to reach the muscle in a significant number of children. Another study presented during the 2013 Annual Meeting of the American Academy of Allergy, Asthma & Immunology (AAAAI) suggested that

2

delivering epinephrine into the muscle allows for more rapid absorption and leads to higher blood levels than if injected into the overlying fat. Considering the rising obesity rates in children, there is concern that epinephrine auto-injectors will not adequately deliver the medication in overweight children who may be experiencing anaphylaxis. Without proper treatment, anaphylaxis can be fatal, therefore it is critical that epinephrine is administered quickly and effectively.

Additional studies have found that patient groups other than children are at risk of not receiving an adequate dose of epinephrine in an anaphylactic emergency. One study published in the American Journal of Emergency Medicine found that the current epinephrine auto injector needle length is inadequate for intramuscular injection, especially among women.

One study conducted at the University of Manitoba assessed absorption of epinephrine when injected subcutaneously and intramuscularly in children. The data showed subcutaneous injections reached the mean maximum plasma concentration at a range of 5 to 120 minutes. Patients injected intramuscularly reached a mean maximum plasma concentration in 8+/-2 minutes. This study supports the conclusion that, in children, recommendations for subcutaneous epinephrine injection are based on anecdotal experience, and should be reevaluated in view of the finding of delayed epinephrine absorption when this route is used. This delay may have important clinical implications during an episode of systemic anaphylaxis. Thus, the intramuscular route of injection is preferable. According to the National Institute of Allergy and Infectious Disease, if epinephrine is not given promptly, rapid decline and death could occur within 30 to 60 minutes. These findings clearly illustrate the need for patients experiencing anaphylactic shock to receive rapid intramuscular injections of epinephrine.

A third disadvantage of auto injectors is the high cost, which particularly affects developing nations.

While vial and syringe combinations for epinephrine administration have the advantage of being economical, they are cumbersome to carry on a daily basis, time-consuming to load and administer, and most importantly, open the door to dosage errors. Patients requiring self-administration are typically not capable of loading a syringe and administering the correct dosage during an anaphylactic event.

Therefore, what is needed is an economical alternative for administering injections during an anaphylactic event that overcomes the disadvantages of prior art injectors.

**BRIEF SUMMARY OF THE INVENTION**

It is an object of the invention to provide an injection device configured to deliver an accurate, personalized, rapid injection of a critical care life-saving dosage of epinephrine when speed, dosage accuracy and proper intramuscular delivery during an anaphylactic event are critical, particularly in children and other small body mass patients.

It is another object of the invention to provide an injection device calibrated to a patient's weight and/or Body Mass Index (BMI) in order to administer a precise dosage with the most operative needle length to improve the speed, accuracy, safety and effectiveness of the administration of a critical care injection such as epinephrine, or other medication.

It is another object of the invention to provide an injection device that is cost-effective and easy to use.

It is another object of the invention to provide an injection device configured to administer an injection in one continu-

3

ous motion using the energy from the force of injection, and not stored energy within the device.

It is another object of the invention to provide an injection device that operates by way of sequential "stages" in one continuous motion to administer a medical injection from the force of applied pressure against the injection site.

It is another object of the invention to provide an injection device in which a predetermined amount of pressure against the injection site initiates the device to, in sequence, puncture the injection site, penetrate hypodermically to a patient specific depth, and deliver a patient specific dosage from within the reservoir of the device.

It is another object of the invention to provide an injection device configured to deliver an injection by way of momentum initiated by the retraction or other actuation of the needle guard or like component of the device applied by force against the injection site.

It is another object of the invention to provide an injection device for administering a medicinal injection that operates by utilizing momentum from the force of injection to deliver the medication.

To achieve the foregoing and other objects and advantages, in a first embodiment the present invention provides a force actuated injection device including a body adapted to contain medication, a grip slidably disposed over the body, a lock disposed between the body and the grip, the lock when locked preventing movement of the grip relative to the body and when unlocked permitting movement of the grip relative to the body, a needle guard retractable relative to the body, a biasing member biasing the needle guard away from the body, and a resistance band disposed between the body and the needle guard adapted to resist initial needle guard retraction relative to the body, wherein force on the needle guard exceeding a predetermined amount overcomes a frictional force of the resistance band allowing the needle guard to retract relative to the body to unlock the grip relative to the body such that movement of the grip relative to the body causes medication to be delivered from the body.

In one aspect, the predetermined amount of force required to overcome the frictional force of the resistance band may be greater than an amount of force required to overcome the biasing member.

In another aspect, the resistance band may be positioned surrounding one end of the body between an outer surface of the body and an inner surface of the needle guard and may frictionally engage the inner surface of the needle guard.

In another aspect, the resistance band may resist only initial retraction of the needle guard relative to the barrel.

In another aspect, the needle guard may include a tubular body having at least one axially extending member that travels along a surface of the body to engage the lock.

In another aspect, the at least one axially extending member may drive the lock radially inward out of engagement with the grip to unlock the grip relative to the body.

In another aspect, the device may include a pair of diametrically opposed locks and the needle guard comprising a pair of diametrically opposed axially extending members each adapted to engage one of the pair of diametrically opposed locks.

In another aspect, the force actuated injection device does not store injection energy.

In another aspect, a window may be formed through the grip.

In another aspect, the body may include diametrically opposed guide rails for guiding diametrically opposing members extending from one end of the needle guard.

4

In another aspect, the lock may be biased radially outward toward the grip and engages the grip when locked.

In another aspect, the body includes an internal fluid reservoir.

In another aspect, the body is adapted to hold a syringe and the grip is adapted to engage a syringe plunger.

Embodiments of the invention can include one or more or any combination of the above features, aspects and configurations.

Additional features, aspects and advantages of the invention will be set forth in the detailed description which follows, and in part will be readily apparent to those skilled in the art from that description or recognized by practicing the invention as described herein. It is to be understood that both the foregoing general description and the following detailed description present various embodiments of the invention, and are intended to provide an overview or framework for understanding the nature and character of the invention as claimed. The accompanying drawings are included to provide a further understanding of the invention and are incorporated in and constitute a part of this specification.

#### BRIEF DESCRIPTION OF THE DRAWINGS

Features, aspects and advantages of the present invention are better understood when the following detailed description of the invention is read with reference to the accompanying drawings, in which:

FIG. 1 is an isometric view of an injection device according to a first embodiment of the invention;

FIG. 2 is an exploded view of the device of FIG. 1;

FIG. 3 is an isometric view of the barrel of the device of FIG. 1;

FIG. 4 is a longitudinal sectional view through the barrel of FIG. 3;

FIG. 5 is a cutaway view of the barrel of FIG. 3;

FIG. 6 is an isometric view of the plunger assembly of the device of FIG. 1;

FIG. 7 is a cutaway view of the plunger assembly of FIG. 6;

FIG. 8 is an exploded view of the plunger assembly of FIG. 6;

FIG. 9 is an isometric view of the needle guard of the device of FIG. 1;

FIG. 10 is an isometric view of the plunger stop of the device of FIG. 1;

FIG. 11 is an isometric view of the plunger lock of the device of FIG. 1;

FIG. 12 is an isometric view of the removable safety cap of the device of FIG. 1,

FIG. 13 is an isometric view of the device of FIG. 1 shown with the plunger lock in the non-operative position;

FIG. 14 is an isometric view of the device of FIG. 13 shown from another angle;

FIG. 15 shows a medication vial aligned with the needle guard;

FIG. 16 shows the needle inserted into the vial;

FIG. 17 shows the device of FIG. 1 loaded with a predetermined dosage of medication and with the plunger lock in the operative position;

FIG. 18 is a detailed view of the plunger lock;

FIG. 19 is a detailed view of the plunger lock and plunger stop engagement;

FIG. 20 shows the device of FIG. 1 with the needle guard retracted and with the plunger in the starting position;



5

FIG. 21 is a detailed view illustrating movement of the needle guard triggering release of the plunger;

FIG. 22 shows the device of FIG. 1 with the needle guard retracted and the plunger fully depressed;

FIG. 23 is an isometric view of an injection device according to a second embodiment of the invention;

FIG. 24 is a longitudinal sectional view through the device of FIG. 23;

FIG. 25 shows the device of FIG. 23 with the safety cap removed;

FIG. 26 shows the device of FIG. 23 with the needle guard partially retracted;

FIG. 27 is an exploded view of the device of FIG. 23;

FIG. 28 is another exploded view of the device of FIG. 23;

FIG. 29 is a sectional view through the device of FIG. 23 showing the device loaded and ready for use;

FIG. 30 is a sectional view through the device of FIG. 23 showing the needle guard partially retracted and the plunger locked;

FIG. 31 is a sectional view through the device of FIG. 23 showing the needle guard fully retracted and the plunger released;

FIG. 32 is an isometric view of an injection device according to a third embodiment of the invention;

FIG. 33 shows the device of FIG. 32 with the safety cap removed and the plunger lock in the non-operative position;

FIG. 34 shows the vial aligned with the needle guard;

FIG. 35 shows the needle inserted into the vial;

FIG. 36 shows the device loaded with a predetermined dosage;

FIG. 37 shows the loaded device with the plunger lock in the operative position;

FIG. 38 is a detailed view of the plunger lock;

FIG. 39 shows the needle depth setting of the device;

FIG. 40 shows the needle depth adjuster;

FIG. 41 is a detailed view of the needle depth adjuster;

FIG. 42 shows the device of FIG. 32 loaded and set to a predetermined needle depth;

FIG. 43 shows the needle guard retracted and the plunger in the starting position;

FIG. 44 shows the needle guard retracted and the plunger fully depressed;

FIG. 45 is a sectional view through the device showing plunger release;

FIG. 46 is an exploded view of the device of FIG. 32;

FIG. 47 is an isometric view of the barrel of the device of FIG. 32;

FIG. 48 is a longitudinal sectional view through the barrel;

FIG. 49 is a cutaway view of the barrel of FIG. 47;

FIG. 50 is an isometric view of the plunger assembly of the device of FIG. 32;

FIG. 51 is a cutaway view through the plunger assembly;

FIG. 52 is an exploded view of the plunger assembly of FIG. 50;

FIG. 53 is an isometric view of the needle guard of the device of FIG. 32;

FIG. 54 is an isometric view of the plunger stop of the device of FIG. 32;

FIG. 55 is another isometric view of the plunger stop;

FIG. 56 is an isometric view of the plunger lock of the device of FIG. 32;

FIG. 57 is an isometric view of an injection device according to a fourth embodiment of the invention;

FIG. 58 is a longitudinal sectional view through the device of FIG. 57;

6

FIG. 59 shows the device of FIG. 57 loaded with a predetermined dosage and ready for use;

FIG. 60 shows the device of FIG. 57 with the needle guard partially retracted;

FIG. 61 is an exploded view of the device of FIG. 57;

FIG. 62 is another exploded view of the device of FIG. 57;

FIG. 63 is a sectional view through the device of FIG. 57 showing the needle guard extended and the plunger assembly locked;

FIG. 64 is a sectional view through the device of FIG. 57 showing the needle guard partially retracted and the plunger assembly locked;

FIG. 65 is a sectional view through the device of FIG. 57 showing the plunger assembly unlocked;

FIG. 66 is an isometric view of an injection device configured to hold a pre-filled syringe;

FIG. 67 is a vertical cross-section through the injection device of FIG. 65;

FIG. 68 shows the pre-filled syringe removed from the injection device;

FIG. 69 shows a syringe lock in an unlocked position;

FIG. 70 shows the syringe lock in a locked position;

FIG. 71 shows the needle guard extended thereby concealing the needle;

FIG. 72 shows the needle guard retracted and the plunger fully depressed;

FIG. 73 is an exploded view of the injection device of FIG. 66; and

FIG. 74 is another exploded view of the injection device of FIG. 66.

#### DETAILED DESCRIPTION OF THE INVENTION

The present invention will now be described more fully hereinafter with reference to the accompanying drawings in which exemplary embodiments of the invention are shown. However, the invention may be embodied in many different forms and should not be construed as limited to the representative embodiments set forth herein. The exemplary embodiments are provided so that this disclosure will be both thorough and complete, and will fully convey the scope of the invention and enable one of ordinary skill in the art to make, use and practice the invention. Like reference numbers refer to like elements throughout the various drawings.

Referring to FIGS. 1-22, a first embodiment of an injection device for administering a dosage of epinephrine or other critical care injection, referred to generically herein as the "medication," is shown generally at reference numeral 100. Although the injection device shown and described herein utilizes a needle for delivering the medication, it is envisioned and intended that the features and configurations described herein may be adapted to and implemented with other types of injection devices for delivering medication hypodermically.

The injection device 100 generally includes a barrel 102, a plunger assembly 106, and a removable safety cap 158. The barrel defines an internal reservoir for containing a volume of medication, also referred to herein as the "dosage." The plunger assembly 106 is axially movable relative to the barrel 102 in a first direction to load the reservoir, and in the opposing direction to force the medication from the reservoir.

The injection device further includes a hypodermic needle 108 configured to attach to one end of the barrel, for example by way of a Luer-lock, to deliver the medication intramuscularly, subcutaneously, or another route, and a retractable

needle guard **110** for guarding the needle and triggering the delivery mechanism of the injection device, as described below. The plunger assembly **106** generally includes a grip **128** that axially overlaps a portion of the barrel **102**, and a plunger **112** configured to be fixed relative to the grip **128**. The plunger **112** may be received through one end of the grip **128** by way of a removable cap **130**. The cap **130** secures in place over one end of the plunger **112** to fix the inner plunger relative to the outer grip. In use, the grip **128** and plunger **112** move together axially relative to the barrel **102** to load and administer the medication. A plunger stop **142** is arranged to control axial movement of the plunger **112** relative to the barrel **102**, and a plunger lock **150** functions to lock the axial position of the plunger assembly **106** relative to the barrel **102**, for example, to set the dosage.

The barrel **102** is an elongate, generally cylindrical body defining an internal passageway therethrough. The hypodermic needle **108** attaches to one end of the barrel **102**, and the opposing end of the barrel is open to receive the internal plunger **112** of the plunger assembly **106**. An airtight seal is provided between one end of the plunger **112** and the inner wall of the barrel reservoir **104** such that the medication is forced from the reservoir through the hypodermic needle as the plunger is axially advanced into the barrel **102**, and cannot pass between the plunger **112** and the reservoir wall. The barrel **102** has a generally cylindrical cross-section with diametrically opposed flat sides. The flat sides may serve to locate indicia **121**, for example, volumetric indicia. Volumetric indicia may include, but is not limited to, mL, IU, CC, MG, MCG, etc. Indicia may also be provided in the form of KG, LB, Body Surface Area (BSA), Body Mass Index (BMI), etc.

The barrel **102** defines a lateral flange **114** adjacent the needle end to facilitate gripping and manipulating the barrel to introduce air into the reservoir. The barrel **102** defines an elongate annular passage **116** for receiving the retracting needle guard **110**. The barrel **102** defines diametrically opposed guide channels **118** for aligning and guiding the needle guard **110** as the needle guard retracts into the barrel. An annular chamber **120** surrounds the end of the reservoir **104** opposite the needle end and receives the plunger stop **142**. The end of the reservoir opposite the needle end terminates in an internal cylindrical stem **122** having diametrically opposed external rails **124** and optional external threading **126** around at least a portion of the circumference thereof, the rails and optional threading functioning, at least in part, to guide the plunger stop into place **142** and facilitate alignment of the plunger stop with the retracting needle guard **110**.

The plunger assembly **106** generally includes the outer grip **128** and the internal plunger **112**. The grip **128** has a cross-section similar to that of the barrel **102**, albeit a larger diameter, such that a portion of the grip axially overlaps the barrel and rotation therebetween is prevented. The plunger **112** terminates at one end in a seal **117** and at the opposing end in an enlarged head **118** trapped between an annular flange of the grip **128** and the removable cap **130**. The grip **128** defines one or more windows, with window **132** indicating the dosage and defined along a flat face of the grip corresponding in alignment with the underlying volumetric indicia on the barrel. Windows **134** are defined thru an arcuate portion of the grip **128** and are used to accessing and manipulate the plunger lock therethrough. The grip **128** may define exterior gripping features **136** to facilitate gripping. An enlarged lateral flange **138** functions as a palm rest during the injecting motion.

The needle guard **110** is a generally cylindrical body having diametrically opposed guide rails **140** extending axially from one end thereof. The needle guard **110** and plunger assembly are calibrated to resist a specific amount of force to actuate the needle guard “upwards” (i.e., in the direction of the barrel). The needle guard **110** “unlocks” the delivery mechanism upon the applied pressure from the injection site, which in turn allows the medication to be delivered. The needle guard **110** retracts within the barrel as the hypodermic needle advances into the patient. At the same time, the guide rails **140** are guided along the guide channels within the barrel. The needle guard **110** retracts as it comes into contact with the injection site from the force of the injection. The needle guard **110** may have a predetermined length thereby providing a predetermined amount of linear length between the needle guard end and the needle tip, thus allowing a calibrated amount of momentum to be in place for the needle tip to penetrate the patient. Prior to injection, the needle guard **110** also functions to conceal the needle and prevent accidental sticks. The needle guard **110** also desensitizes the injection site just prior to the needle puncturing the skin lessening/reducing the initial pricking/pinching sensation felt by the patient. A relatively tight fit may be provided between the barrel **102** and the needle guard **110** such that a predetermined force is needed to move the needle guard relative to the barrel. With such a tight fit, the needle guard is not free to move simply by changing the orientation of the tool, but requires that the needle guard either come into contact with the injection site or is otherwise purposefully handled and moved (e.g., retracted or pulled).

The plunger stop **142**, which sits within the annular chamber around the stem, generally includes an annular band **144** having diametrically opposed arms **146** that extend in the same direction away and outward from the band. The band defines inner channels **148** that align with the stem rails to align the plunger stop **142** relative to/on the stem. The arms **146** are resiliently deformable and biased radially outward such that when the needle guard **110** is not fully retracted the ends of the arms engage within a slot in the plunger lock to prevent axial movement of the plunger relative to the barrel, and are forced radially inward and out of engagement with the plunger lock when the needle guard is retracted, thereby allowing the plunger to move axially relative to the barrel to deliver the medication.

The plunger lock **150** functions to lock the axial position of the plunger assembly relative to the barrel, for example, to set the dosage. The plunger lock **150** includes diametrically opposed finger grips **152** that are accessible through the windows located along the sides of the grip. The plunger lock **150** may be made from a resiliently deformable material capable of being compressed radially inward in order to disengage teeth **154** from corresponding teeth located on the inner wall of the grip (see FIG. **18** at **156**). The plunger lock **150** installs around the stem of the barrel, and the arms of the plunger stop rest against the end or a shoulder of the plunger lock **150** when the plunger lock is in the operative position, thus preventing axial movement between the plunger and barrel. As the arms **146** are urged radially inward in response to force from the retracting needle guard **110**, the arms move out of engagement with and clear of the plunger lock **150**, thereby allowing the plunger assembly to move axially in the direction of the barrel to deliver the medication.

The needle cap **158** is removed before syringe use. The needle cap **158** defines an outer cylindrical portion **160** that surrounds at least a portion of the needle guard, and an inner cylindrical needle holder **162** that surrounds the hypodermic

9

needle. The needle holder **162** may function to hold the hypodermic needle during coupling with the Luer-lok of the barrel.

FIGS. **13-22** detail the sequential steps of readying, loading and administering an injection utilizing injection device **100**. First, referring to FIG. **13**, with the safety cap removed, the plunger lock **150** is moved to a predetermined starting position at the end of the grip **128** opposite the barrel, thereby allowing the plunger assembly to move freely relative to the barrel. Next, referring to FIG. **14**, the user grips the grip **128**, resting the top of the hand firmly against the underside of the grip flange **138** and placing the thumb against the underside of the barrel flange in order to push the barrel **102** in the direction away from the plunger to introduce air into the reservoir. Next, referring to FIG. **15**, the user aligns the vial with the needle guard flange. Next, referring to FIG. **16**, the user pushes the injection device or vial until the needle penetrates the vial and the “-0-” or equivalent zero dosage shows through the dosage window **132**. The plunger is then withdrawn from the barrel **102** to fill the reservoir to the desired dosage amount that appears through the dosage window **132**.

Next, referring to FIG. **17**, with the vial removed, the needle guard **110** is pulled out to a starting position and the plunger lock is “pinched” to compress the plunger lock radially to slide it into the operative “locked” position against the arms of the plunger stop. FIG. **18** details the plunger lock **150** engaging the teeth **156** on the inner surface of the grip. FIG. **19** details the plunger lock **150** engagement with the barrel and with the arms of the plunger stop.

Next, referring to FIG. **20**, to administer the injection, the user grips the grip **128** seating the hand against the grip flange **138**. The user applies a downward pressure on the grip and grip flange **138**, thereby causing force between the injection site and needle guard **110** until the needle guard collapses and retracts upward into the barrel **102** a predetermined distance and the needle penetrates the injection site in response to forward pressure against the injection site. Referring to FIG. **21**, as the needle guard **110** retracts, the needle guard rails **140** move into contact with the plunger stop arms **146**. As the needle guard **110** continues to retract, the contact between the rails **140** and the arms **146** urges the arms radially inward, thereby releasing the plunger stop and allowing the plunger to move to complete the injection. Referring to FIG. **22**, as part of the continuous injecting motion, the plunger continues its “downward” motion to force the medication from the reservoir and out through the needle **108**. Thus, needle guard retraction, needle insertion, plunger release, and plunger movement are sequentially performed in one continuous motion from force against the injection site, and plunger movement is stayed until the needle is fully inserted into the patient to a predetermined depth. After the injection is complete, the injection device may be discarded or reused after sterilizing the device and replacing the needle.

Referring to FIGS. **23-31**, a second embodiment of an injection device for administering a dosage of medication is shown generally at reference numeral **200**. Like injection device **100**, injection device **200** generally operates by way of user force and momentum to sequentially insert the needle, trigger release of the plunger, and deliver the medication, in one continuous motion. Unlike injection device **100**, in which the dosage can be adjusted, injection device **200** is pre-loaded with a predetermined dosage. The dosage can be customized based on weight, BMI and BSA, among other factors.

10

The injection device **200** includes a needle cap **258** removed prior to use. The needle cap **258** includes an inner cylindrical needle holder **259** that conceals the hypodermic needle **108**. The needle cap **258** may be held in place, for example, by way of interference fit around a portion of the grip **228**. Injection device **200** generally includes a barrel **202** defining an internal reservoir **204** containing a predetermined volume of medication, a plunger assembly **206** axially movable relative to the barrel to force the medication from the reservoir, a hypodermic needle **108** for delivering the medication intramuscularly, subcutaneously, or other route, and a retractable needle guard **210** for releasing the plunger assembly **206** to allow axial movement of the plunger relative to the barrel.

The barrel **202** may be a unitary or multi-piece body. The hypodermic needle **108** attaches to one end of the barrel **202**, and the opposing end of the barrel is open to receive an internal plunger **212** of the plunger assembly **206** therein. An airtight seal is provided between the plunger **212** and the reservoir **204** such that the pre-loaded medication is forced through the needle as the plunger moves axially in the direction of the barrel **202**. A portion of the barrel **202** has a generally cylindrical cross-section with diametrically opposed flat sides, which correspond in shape to the overlying grip **228** to prevent rotation therebetween when axially overlapped. The end of the reservoir adjacent the needle, i.e., the interface between the reservoir and needle, may be sealed with a foil or like seal that ruptures in response to the initial plunger movement and/or consequential pressure build-up in the reservoir. The rupture of the seal allows the medication to flow from the reservoir **204** to the needle.

In response to predetermined pressure against the injection site, the needle guard **210** rapidly retracts into the barrel **202**. The barrel **202** defines an annular, internal space **216** for receiving the cylindrical body of the retracting needle guard **210**, and diametrically opposed guide channels **218** for receiving and aligning the corresponding guide rails **240** extending from the needle guard body. The plunger stop **242** positioned around the axially extending stem **222** of the barrel **202** includes diametrically opposed arms **246** extending in the direction of the top of the grip **228**. In use, the arms **246** are resiliently deformable and biased radially outward such that when the needle guard **210** is extended the ends of the arms engage within corresponding slots **260** thru the grip **228** to prevent axial movement of the plunger **212** relative to the barrel **202**. When the needle guard **210** is retracted into the barrel **202** a sufficient distance, the needle guard urges the arms **246** radially inward out of engagement with the grip slots **260** such that the rails **240** can pass beyond the arms allowing the grip **228** to continue its motion in the direction of the injection site, thereby allowing the medication to be delivered. Features **260** facilitate gripping and handling.

Referring to FIGS. **32-56**, a third embodiment of an injection device for administering a dosage of medication is shown generally at reference numeral **300**. Injection device **300** operates much like injection devices **100** and **200**, with the added functionality of adjustable needle length. It should be understood that the adjustable needle length functionality may be added to one or more of injection devices **100** and **200** for the added benefit of patient customization, among other advantages.

Grip **328** includes a dosage window **332** through which the set dosage is indicated, and a needle depth window **334** through which the needle depth is indicated, for example, in units of millimeters. FIG. **32** shows the injection device **100** with the needle cap **358** installed, thereby concealing the

11

underlying needle and preventing contamination. FIG. 33 shows the injection device in the “starting position,” with the dosage set to “-0-”, the plunger lock 350 disengaged allowing the plunger assembly 306 to move freely, and with the minimum needle depth set. Vial or ampoule alignment and plunger lock operation work in the same manner described above with reference to injection device 100.

The needle depth may be adjusted to change the depth of needle penetration. Depth of penetration is controlled by adjusting the point at which the needle guard 310 engages the plunger stop 342 to trigger the release of the plunger 312, and the adjustment is actuated by rotating a dial 370 accessible through the flat face of the grip 328 opposite the face including the dosage and needle depth indication windows. The dial 370 is part of an assembly also including one or more dial nuts 372 threaded onto the plunger stem 322. The dial 370 and one or more dial nuts 372 are meshed gears, wherein rotational motion of the dial 370 causes rotational motion of the one or more dial nuts 372. When turned, the one or more dial nuts 372 force the plunger stop 342 up or down the plunger stem 322 depending on the desired needle depth. In the embodiment shown, the dial assembly includes upper and lower dial nuts 372 positioned vertically above and below the plunger stop 342, and the dial 370 is held in a meshed arrangement with the dial nuts 372 by way of a dial holder 374 positioned adjacent a side of the plunger stop 342. With the dosage and needle depth set, the injection device 300 operates to deliver the medication similar to injection device 100.

The barrel 302 generally includes the reservoir 304 for holding the medication, the stem 322 for locating the plunger stop 342, and the various guides and chambers for locating and guiding the retractable needle guard 310. Teeth on the stem 322 serve to grip the barrel 302 to prevent the plunger from pulling out of the barrel once the medication has been loaded and the plunger lock set to the predetermined dosage amount for the injection. External threads 376 on the stem 322 engage the internal threading on the one or more dial nuts 372 to adjust the needle depth in response to rotational movement of the dial 370. The grip flange 338 serves as a palm rest during the injecting motion.

The plunger stop 342 includes resilient arms 346 for engaging the plunger lock 350 (see FIG. 56) to prevent the plunger 312 from advancing into the reservoir 304 until released. Needle depth indicia may be located on a face of the plunger stop visible through the needle depth window 334. A slot is provided on one face of the dial for attaching the dial holder.

Referring to FIGS. 57-65, a fourth embodiment of an injection device for administering a dosage of epinephrine or other critical care injection is shown generally at reference numeral 400. Injection device 400 generally includes a barrel 402, plunger assembly 404, needle guard 406 and removable safety cap 408. As best shown in FIG. 58, the grip 410 and internal plunger 412 may be integrally formed, and the grip 410 may have a pistol grip shape to facilitate gripping.

The needle guard 406 is biased in the direction concealing the hypodermic needle 414. Optional helical spring 416 may be maintained between the needle guard 406 and the end of the barrel 402, requiring a small amount of force to overcome the spring force to retract the needle guard into the barrel. In an alternative embodiment, a resistance band 422 may be positioned surrounding the end of the barrel between the outer surface of the barrel and inner surface of the needle guard 406. The resistance band resists movement of the needle guard 406 relative to the barrel 402, requiring a

12

predetermined degree of force to initiate retraction movement of the needle guard into the barrel. The helical spring 416 may be optionally included to provide smooth movement of the needle guard 406 relative to the barrel 402 throughout its range of movement. In one particular embodiment, both a resistance band and helical spring may be utilized, with the former resisting initial needle guard movement until a predetermined threshold resistance force is overcome, and the later providing smooth needle guard movement thereafter and biasing force in the direction concealing the needle.

A pair of axially extending spaced arms 418 of the needle guard 406 travel along the inner surface of the barrel 402. As the needle guard 406 is retracted into the barrel 402, the ends of the arms 418 engage a pair of spaced locks 420 to drive them radial inward to unlock the plunger assembly 404 from the barrel, thereby allowing the plunger 412 to be axially advanced into the reservoir.

The spaced locks 420 are arranged between the barrel 402 and inner surface of the grip 410 and are spring biased in the direction of the grip. Each lock 420 is ramped such that continued movement of the arms 418 in the direction of the locks causes the locks to be driven radially inward out of contact with the grip 410. In the “locked” position, each lock sits within a recess defined along the inner surface of the grip 410. In the “unlocked” position, the locks are driven from their respective recesses such that the plunger assembly 404 is free to move relative to the barrel 402. FIG. 63 shows the needle guard 406 fully extended and the locks 420 in their locked position. FIG. 64 shows the needle guard 406 partially retracted into the barrel and with the arms 418 making initial contact with the locks 420. FIG. 65 shows the needle guard 406 retracted a sufficient amount such that the locks 420 are driven urged radially inward and out of contact with the grip 410.

Referring to FIGS. 66-74, a device for controlling actuation of a pre-filled syringe is shown generally at reference numeral 500. The device 500 generally operates to hold a syringe filled with a predetermined medication, i.e., a pre-filled syringe, and actuate the syringe to administer the predetermined medication in response to pressure against the injection site as discussed above in connection with the foregoing embodiments. In other words, the injection device 500 is momentum driven in that force against the injection site causes the needle guard of the device to partially retract into the barrel holder a predetermined distance to unlock the grip from the barrel holder to allow the grip to move in a direction of the injection site to advance the syringe plunger into the syringe barrel to deliver a dosage of the predetermined medication in one continuous motion.

The injection device 500 generally includes a barrel holder 502 having an internal axial passageway adapted to receive/hold a barrel 504 of the pre-filled syringe 506, a grip 508 slidably disposed over the barrel holder and adapted to engage a plunger 510 of the pre-filled syringe, and a needle guard 512 configured to partially retract into the barrel holder in response to force against the injection site to unlock the grip from the barrel holder. When locked, the grip 508 is fixed relative to the barrel holder 502 thereby preventing axial movement therebetween and consequently plunger movement relative to the syringe barrel. When unlocked, the grip 508 is able to move relative to the barrel holder 502 thereby pushing the plunger 510 in the direction of the barrel 504 to drive the predetermined medication out through the needle 514. FIG. 66 shows the injection device 500 locked and with a removable cover 516 installed, while FIG. 67 is a vertical sectional view through the injection

13

device shown with the cover removed and with the grip **508** locked relative to the barrel holder **502**.

The grip **508** has a removable cap **518** at one end thereof adapted to engage an end of the plunger **510** when installed on the grip. The removable cap **518** has an opening **520** through a sidewall **522** thereof for allowing the removable cap to slide over or past the end of the plunger **510** when installing the removable cap on the grip **508**. As best shown in FIG. **68**, the pre-filled syringe **506** is installed into the injection device **500** by first removing the cap **518** from the grip **508**, next inserting the syringe **506** into the injection device, and finally replacing the cap **518** by sliding the cap onto the grip **508** from the side of the grip. The cap **518** and grip **508** can secure together through slide-lock engagement or equivalent manner of engagement. When the cap **518** is installed on the grip **508** the inner face of the end of the cap can directly contact the end of the enlarged head of the plunger.

As best shown in FIGS. **69** and **70**, the syringe barrel locks to the barrel holder **502** by way of a lock **524** positioned at one end of the barrel holder configured to lock the syringe barrel relative to the barrel holder. As shown, the finger flange **526** formed at the end of the syringe barrel sits within a recess formed in the end of the barrel holder **502**. The lock **524**, which slides along a rail **528** laterally adjacent the recess, defines a feature at one end that engages through a surrounding relationship one wing of the finger flange. More specifically, the lock **524** captures the wing of the finger flange between two arms thereof such that sliding the lock **524** along the rail **528** causes the finger flange to turn within the recess. The recess can be topped with a flange that allows the finger flange to pass thereunder when rotated to the locked position, thereby compressing the finger flange and locking the syringe barrel relative to the barrel holder. In the locked position, the syringe barrel is prevented from moving axially or rotationally relative to the barrel holder. The lock **524** is unlocked by sliding the lock in the opposite direction along the rail **528**.

As best shown in FIG. **66**, the injection device **500** can further include a window **530** through a sidewall thereof for viewing indicia **532** or fluid of the pre-filled syringe installed within the device to ensure proper installation among other purposes.

Referring to FIGS. **71** and **72**, the needle guard **512** partially retracts into the barrel holder in response to force against the injection site. FIG. **71** shows the needle guard **512** fully extended thereby concealing the needle. FIG. **72** shows the needle guard **512** partially retracted thereby exposing the needle. As can be seen referring to FIGS. **67**, **73** and **74**, grip **508** locking relative to the barrel holder **502** is achieved using a pair of diametrically opposed grip locks **534** outwardly biased in the direction away from the barrel holder **502**. Biasing can be achieved using helical springs **536** or equivalent biasing member positioned inward of the grip locks **534**. Each grip lock **534** seats within a dedicated opening or recess **538** defined in the sidewall of the barrel holder **502**. When assembled, the grip locks **534** are biased outward to engage an inner surface of the grip **508** to lock the grip to the barrel holder **504**, and are arranged to move out of contact with the grip **508** in response to contact with the needle guard **512** as the needle guard retracts into the barrel holder **502**.

The needle guard **512** generally includes a cylindrical tubular body **540** having a pair of diametrically opposed arms **542** extending axially from one end of the cylindrical body in a direction of the grip **508**. As the needle guard **512**, and particularly the arms **542**, retract into the barrel holder

14

**502** a predetermined distance, the arms **542** contact the grip locks **534** thereby urging the grip locks inward and out of contact with the grip **508** to unlock the grip from the barrel holder **502** to permit axial movement therebetween. The arms **542** can be guided along elongate channels **544** within the interior of the barrel holder **502**. A biasing member, such as a helical spring **546** can be disposed between one end of the barrel holder **502** and the needle guard **512** arranged to bias the needle guard in the direction away from the barrel holder, thereby requiring force against the injection site sufficient to overcome the spring force in order to allow the needle guard **512** to retract.

As in the above injection device embodiments, the injection device **500** can further include a needle depth adjustment mechanism including a plunger stop, a rotating dial accessible through a face of the grip, and one or more dial nuts adapted to thread on the plunger, the rotating dial adapted to mesh with the one or more dial nuts such that rotation of the dial rotates the one or more dial nuts to move the plunger stop up or down the plunger.

In any of the foregoing embodiments, the injection device may be included in a kit for increasing the effectiveness and safety of an injection. The kit may include the injection device marked with patient weight indicia, for example from 1-40 Kg for pediatric patients, or other patient characteristic indicia. The kit may further include one or more of a BMI card to assess a patient's BMI, correlating needles having a predetermined length to be used with the injection device, and a vial or other medicament container including a predetermined volume to support predetermined dosages.

The injection device or kit may be volumetrically predetermined for the medication and targeted patient class and/or patient segment dosage range. It is contemplated herein that the injection device may be filled by a user to a specific dosage, weight or other measurement, and the dosage stored and "locked" within the reservoir of the injection device for later use. The injection device may be pre-filled with the medicament to a predetermined volume to fulfill the administration requirements of a specific patient class, weight or other feature, as an example, pre-filled for a pediatric patient class of 40 Kg or under, or pre-filled for an infant weight class from 1-10 Kg. The injection device may be pre-filled with medication corresponding to a specific patient weight class or other patient characteristic dosage or dosage range, and predetermined needle length. The injection device may be pre-filled with a predetermined volume of medication to fulfill the administration requirements of a specific patient class, weight, BMI or other patient characteristic, and may also include a patient specific needle or other hypodermic delivery apparatus with packaging, labeling, color or text that would assist a user in a quick identification of the most optimal injection device for a particular patient.

The foregoing description provides embodiments of the invention by way of example only. It is envisioned that other embodiments may perform similar functions and/or achieve similar results. Any and all such equivalent embodiments and examples are within the scope of this disclosure.

What is claimed is:

1. A force actuated injection device, comprising:
  - a body adapted to contain medication;
  - a grip slidably disposed over the body;
  - a lock disposed between the body and the grip, the lock when locked preventing movement of the grip relative to the body and when unlocked permitting movement of the grip relative to the body;
  - a needle guard retractable relative to the body;

15

a biasing member biasing the needle guard away from the body; and

a resistance band disposed between the body and the needle guard adapted to resist initial needle guard retraction relative to the body;

wherein force on the needle guard exceeding a predetermined amount overcomes a frictional force of the resistance band allowing the needle guard to retract relative to the body to unlock the grip relative to the body such that movement of the grip relative to the body causes medication to be delivered from the body; and

wherein the resistance band is positioned surrounding one end of the body between an outer surface of the body and an inner surface of the needle guard and frictionally engages the inner surface of the needle guard.

2. The force actuated injection device according to claim 1, wherein the predetermined amount of force required to overcome the frictional force of the resistance band is greater than an amount of force required to overcome the biasing member.

3. The force actuated injection device according to claim 1, wherein the resistance band resists only initial retraction of the needle guard relative to the body.

4. The force actuated injection device according to claim 1, wherein the needle guard comprises a tubular body having at least one axially extending member that travels along a surface of the body to engage the lock.

5. The force actuated injection device according to claim 4, wherein the at least one axially extending member drives the lock radially inward out of engagement with the grip to unlock the grip relative to the body.

6. The force actuated injection device according to claim 1, the lock further comprising a pair of diametrically opposed locks and the needle guard comprising a pair of diametrically opposed axially extending members each adapted to engage one of the pair of diametrically opposed locks.

7. The force actuated injection device according to claim 1, wherein the force actuated injection device has no stored internal energy for administering an injection.

8. The force actuated injection device according to claim 1, further comprising a window formed through the grip.

9. The force actuated injection device according to claim 1, wherein the body comprises diametrically opposed guide rails for guiding diametrically opposing members extending from one end of the needle guard.

16

10. The force actuated injection device according to claim 1, wherein the lock is biased radially outward toward the grip and engages the grip when locked.

11. The force actuated injection device according to claim 1, wherein the body includes an internal fluid reservoir.

12. The force actuated injection device according to claim 1, wherein the body is adapted to hold a syringe and the grip is adapted to engage a syringe plunger.

13. A force actuated injection device, comprising:

a body adapted to contain medication;

a grip slidable relative to the body;

a lock that when locked prevents movement of the grip relative to the body and when unlocked permits movement of the grip relative to the body;

a needle guard retractable relative to the body;

a biasing member biasing the needle guard away from the body; and

a resistance band disposed between the body and the needle guard frictionally engaging the needle guard;

wherein force on the needle guard exceeding a predetermined amount overcomes a frictional force of the resistance band allowing the needle guard to retract relative to the body to unlock the grip relative to the body such that movement of the grip relative to the body causes medication to be delivered from the body; and

wherein the resistance band is positioned surrounding one end of the body between an outer surface of the body and an inner surface of the needle guard and frictionally engages the inner surface of the needle guard.

14. The force actuated injection device according to claim 13, wherein the predetermined amount of force required to overcome the frictional force of the resistance band is greater than an amount of force required to overcome the biasing member.

15. The force actuated injection device according to claim 13, wherein the resistance band resists only initial retraction of the needle guard relative to the body.

16. The force actuated injection device according to claim 13, wherein the force actuated injection device has no stored internal energy for administering an injection.

17. The force actuated injection device according to claim 13, wherein the body includes an internal fluid reservoir.

18. The force actuated injection device according to claim 13, wherein the body is adapted to hold a syringe and the grip is adapted to engage a syringe plunger.

\* \* \* \* \*